medivon.


User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

One massager. Endless ways to improve your life.

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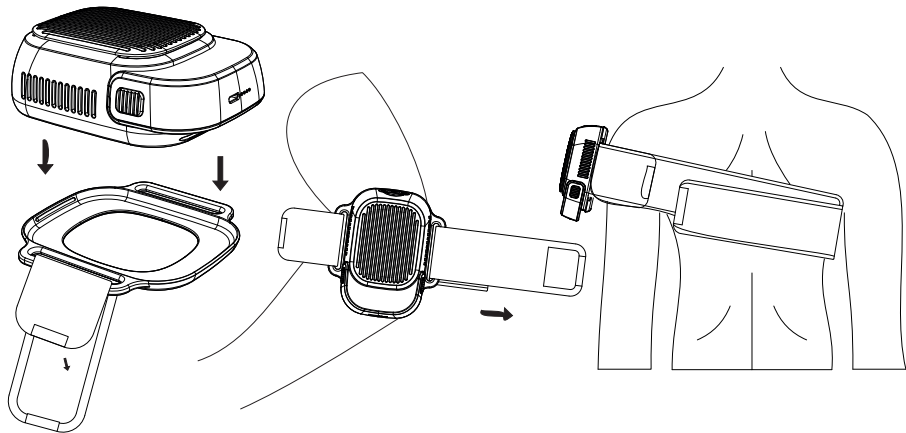The device is designed for localized application of heat, cold, or contrast therapy (alternating cold/warmth) to relieve muscle pain, joint stiffness, and aid recovery after exercise. It can be used directly in the hand or attached to a specific body part with straps.

To begin, hold the cold (snowflake) button for cooling therapy, or the heat (heat wave) button for heat therapy. To activate contrast mode, hold the center contrast button.

Cycle to change the temperature between: cold 16°C / 12°C / 8°C (default 16°C), and heat 35°C / 39°C / 43°C (default 35°C). The therapy runs automatically for 18 minutes (cold), 24 minutes (warm), or 20 minutes (contrast: 5 min cold / 5 min warm x 2), after which the device turns off.

For handheld use, place the device against your body—preferably directly on your skin. For hands-free use, secure the device with the straps: remove it from the frame, thread the straps through the side buckles, wrap them around the desired body part (the logo must be facing out), tighten the Velcro, and adjust the tension. Do not use around your neck; always use both straps. To charge the device, connect the USB-C cable to the USB-A port and charger. While charging, the LED flashes—a blue light indicates a full charge.

Do not use the device while charging. Always turn off the device before cleaning. Use a damp microfiber cloth for cleaning—do not immerse, use detergents, or machine wash. Store in a dry place, away from moisture and sunlight. Do not use the device if you have a pacemaker, implants, heart disease, epilepsy, cancer, infections, pregnancy, open wounds, neurological conditions, or if you are a child or dependent.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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The manufacturer hereby declares that the Medivon DualTherm device complies with, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87, 00-844 Warszawa, Poland.
Made in China.
Marka: Medivon
Model: Medivon DualTherm
email: serwis@medivon.pl
SN:06255904119286445

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>