

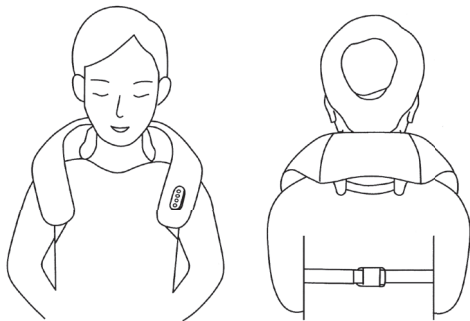
medivon.

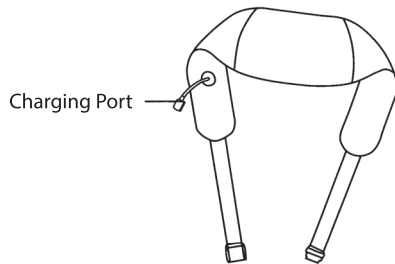User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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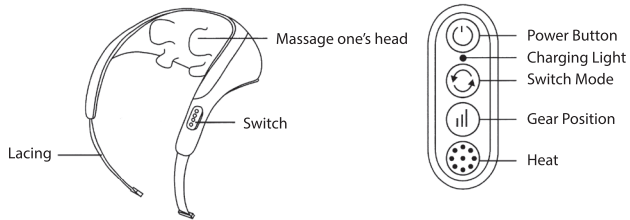Designed for massaging the neck, nape, and upper back, this device combines a kneading massage technique with a heating function. It helps reduce muscle tension, relieves pain caused by prolonged sitting, and improves circulation and flexibility. Its compact, ergonomic design makes it perfect for use at home, at work, or on the go. The device features four massage modes and three intensity levels—changeable by pressing the power button. The heating function activates automatically in each mode. The entire massage cycle lasts approximately 10 minutes, after which the device automatically turns off.

During use, simply rest your neck on the device in a sitting or semi-reclining position. The massage is hands-free—the kneading rollers operate independently. The built-in 1800 mAh battery allows for wireless operation and is charged via a USB-C (5V) port. The LED flashes red while charging and glows green when fully charged. Do not use the device while charging.

The outer protective fabric is removable and can be hand washed in lukewarm water. The casing should be cleaned gently with a damp microfiber cloth. Do not immerse the device in water, machine wash, or use chemicals. Store in a dry place, away from moisture and sunlight.

Do not use the device if you have a pacemaker, cancer, infection, neurological disease, epilepsy, pregnancy, open wounds, or if you are dependent on someone without supervision. The device is not suitable for children. If you experience pain, discomfort, or dizziness, discontinue use immediately. Do not use the massager for more than 10 minutes at a time. If in doubt, consult a doctor.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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The manufacturer hereby declares that the Medivon Volt device complies with, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87, 00-844 Warszawa, Poland.
Made in China.
Marka: Medivon
Model: Medivon Volt
email: serwis@medivon.pl
SN:07255904119286544

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>