nenob[®] beauty

Auri

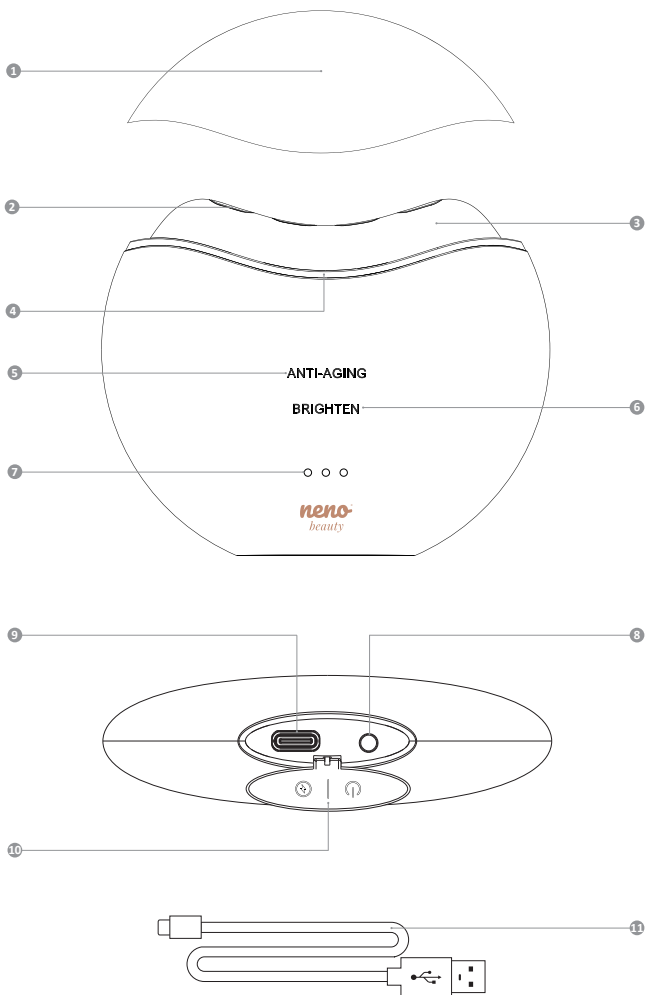
Table of contents

Spis treści

4	instrukcja obsługi	PL
6	User manual	EN
7	Bedienungsanleitung	DE
9	Manuel de l'utilisateur	FR
11	Návod k použití	CZ
13	Používateľská príručka	SK
15	Használati utasítás	HU
16	Användarhandbok	SE
18	Käyttäjän käsikirja	FI
20	Brukerhåndbok	NO
22	Brugermanual	DK
23	Gebruikershandleiding	NL
25	Manual del usuario	ES
27	Manuale utente	IT
29	Manual de utilizare	RO
31	Korisnički priručnik	HR/BA
33	Uputstvo za upotrebu	RS/ME
34	Uporabniški pravilnik	SI
36	Εγχειρίδιο χρήστη	GR
38	Lietotāja rokasgrāmata	LV
40	Vartotojo vadovas	LT
42	Kasutusjuhend	ET
43	Посібник користувача	UA
45	Ръководство за потребителя	BG
47	Упатство за употреба	MK
49	Manual do utilizador	PT

A

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Szanowny Kliencie,

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USER MANUAL

Dear Customer,

Thank you for purchasing the Neno Beauty Auri 5in1 face and neck lifting massager. The device has been designed for comprehensive face and neck care - combining five cutting-edge technologies that promote firming, smoothing and refreshing of the skin and help improve facial contours. This user manual contains all the information necessary for correct use. Please familiarise yourself with it and keep it for future reference..

01. PRECAUTIONS

1. Do not use the device on swollen or inflamed skin, or on skin that is damaged, burnt or immediately after recent surgery.
2. The device is not recommended for people with electronic implants (e.g. pacemaker, insulin pump).
3. The device should not be used in cases of active skin diseases (ulcers, psoriasis, eczema, rashes, recurrent herpes), epidermal damage (scars, wounds), hypersensitivity to electromagnetism, pregnancy and lactation, haemophilia/haemorrhagic disorders, or atopic dermatitis (AD).
4. Do not use the massager if you are taking photosensitising medication (e.g. certain antibiotics, retinoids).
5. The device should not be used by children. Adults with physical, sensory or cognitive limitations should only use the device under adult supervision.
6. If you have any doubts about the suitability of the device for a particular health condition (e.g. skin diseases), consult your doctor.
7. Avoid looking directly at operating LEDs.
8. Do not apply too long on one area of skin.
9. Do not use the device for longer than the recommended time. Use longer than indicated can lead to skin irritation.
10. Do not use the device near water or in high humidity conditions such as a bath, shower or sauna.
11. If you feel discomfort during use, discontinue use immediately.
12. Do not use the device while it is charging. Disconnect it from the power supply when it is fully charged.
13. Keep the power source away from sources of heat and water while charging.
14. Do not use the product if any component is damaged or missing.
15. Do not disassemble or repair the product yourself.
16. Do not use accessories or spare parts other than those supplied.
17. This product is not a toy. Do not allow children to play with the device.
18. Do not store the device and its packaging in a place that is easily accessible to children.
19. Do not store the device near a heat source and do not expose it directly to sunlight. Do not immerse the device in water.

02. CONTENTS OF THE KIT

1. Massager
2. Storage bag
3. Microfibre cloth
4. Type-c cable
5. Pin to open the cap
6. Operating instructions

03. DESCRIPTION OF THE DEVICE SEE FIGURE A

1. Magnetic cover
2. Four massage points
3. Contoured massage head
4. LEDs
5. ANTI-AGING mode indicator
6. BRIGHTEN mode indicator
7. Power indicator
8. On/off/change mode button
9. USB charging port
10. Cap
11. Type-c cable

04. OPERATION OF THE DEVICE

1. The device has two modes:
 - ANTI-AGING (rejuvenating): dedicated to comprehensive regeneration, it combines precise RF waves with sonic massage and an ergonomic Gua Sha-shaped head to deeply stimulate collagen production, smooth wrinkles and firm and shape the contours of the face and neck. The mode has two levels of intensity.
 - BRIGHTEN (brightening): an advanced version of the ANTI-AGING mode, enriched with muscle stimulation using EMS pulses and 620 nm red light therapy. EMS activates mimic muscles, improving skin tone and elasticity, resulting in a visible lifting effect and enhancing facial contours. The red LED light further promotes regeneration, helps reduce the appearance of minor discolouration, reduces grey complexion and restores a healthy, radiant appearance. The mode has two intensity levels.

2. Wash and dry your face and neck thoroughly before use. For optimal results and comfort, we recommend using the device with a conductive gel suitable for RF and EMS technology or with a light, oil-free facial cosmetic that does not impede the smooth gliding of the device over the skin.
3. Select one of the modes:
 - **ANTI-AGING** (rejuvenation): to switch on the device remove the magnetic cover or press and hold button A8 (a beep will sound). Then use the same button to set the ANTI-AGING mode (the mode indicator will illuminate yellow or red and the unit will beep). Yellow illumination indicates a gentle rejuvenation mode, while red illumination is an intense rejuvenation mode. The two intensity levels may differ in the degree of sonic vibration, allowing the device's action to be tailored to individual skin needs.
 - **BRIGHTEN**: to switch on the device, remove the magnetic cover or press and hold button A8 (a beep will sound). Then set the BRIGHTEN mode using the same button (the mode indicator will light up yellow or red and the unit will beep). Yellow backlighting indicates a gentle brightening mode, while red backlighting is an intense brightening mode. The two intensity levels differ in the strength of the EMS pulses and the degree of sonic vibration, allowing the device's action to be precisely tailored to individual skin needs.
4. During the massage, move the device slowly over the face and neck, making upward and outward movements, following the natural direction of lymphatic drainage. The device signals the change of area every minute.
5. The device has a built-in 6-minute timer: six minutes after setting one of the modes, the device will switch off. The time does not add up when changing modes during their operation.
6. To switch off the device before 6 minutes have elapsed, the magnetic cover must be applied or button A8 must be pressed and held down.
7. Once the massage is over, you can move on to the next stages of care.
8. The device has a memory mode - when the massager is switched on again, it will automatically start in the last mode used.
9. Use the massager 2-4 times a week (for 6-12 minutes), choosing the mode that suits your needs.

05. CLEANING AND STORAGE

1. After each use, wipe the massage head with a dry or slightly damp cloth.
2. Let it dry in a well-ventilated place, away from direct sunlight and heat sources, then apply the magnetic cover.
3. Do not immerse the device in water and do not use highly corrosive detergents or alcohol.
4. When not in use, the product should be stored in the bag provided.

06. CHARGING THE DEVICE

1. To start charging your device, open the charging port cap with the supplied pin.
2. During charging, the three power indicator LEDs flash successively from left to right. The device cannot be used during charging.
3. The power indicator lights illuminate continuously when the device is fully charged.
4. When the battery level drops below 10%, the power indicator LEDs flash sequentially from left to right, indicating that the device needs to be charged.
5. If an attempt is made to start the unit with a fully discharged battery, a double beep will sound, the power indicator LEDs will flash twice and the unit will automatically switch off.

07. SPECIFICATION

Power supply: 5V 0.29A

Power: 5W

Battery: Lithium-polymer 500mAh

Timer: 6min

Light wavelengths: red light 620±5 nm

Material: ABS + TPU + PC

Dimensions: 90.4 x 86.4 x 26.4 mm

Weight: 104g

08. WARRANTY CARD

The product comes with a 24-month warranty. Warranty conditions can be found on the website: <https://neno.pl/gwarancja>
 Details, contact and service address can be found at: <https://neno.pl/kontakt>
 Specifications and contents are subject to change without notice. We apologise for any inconvenience.

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