nenob[®] beauty

Luve

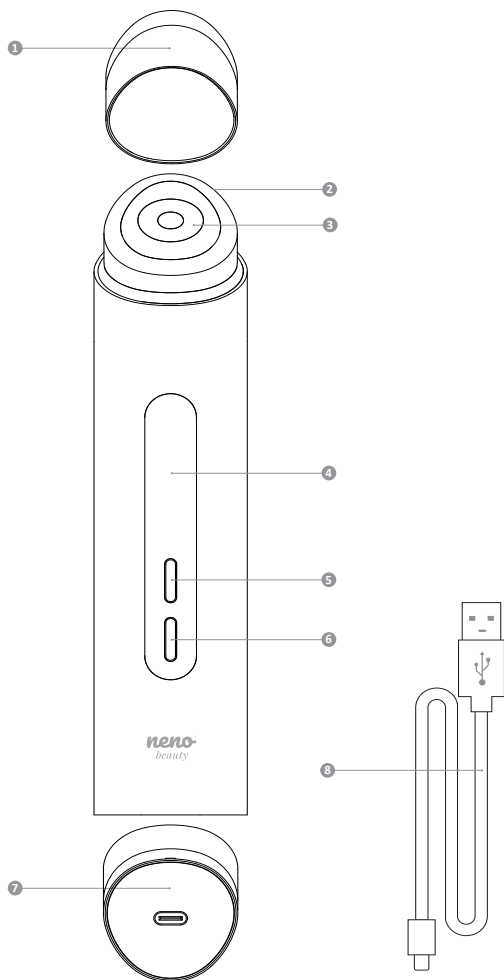
Table of contents

Spis treści

4	Instrukcja obsługi	PL
6	User manual	EN
9	Bedienungsanleitung	DE
12	Manuel de l'utilisateur	FR
14	Návod k použití	CZ
17	Používateľská príručka	SK
20	Használati utasítás	HU
22	Användarhandbok	SE
25	Käyttäjän käsikirja	FI
27	Brukerhåndbok	NO
30	Brugermanual	DK
32	Gebruikershandleiding	NL
35	Manual del usuario	ES
38	Manuale utente	IT
40	Manual de utilizare	RO
43	Korisnički priručnik	HR/BA
45	Uputstvo za upotrebu	RS/ME
48	Uporabniški pravilnik	SI
51	Εγχειρίδιο χρήστη	GR
54	Lietotāja rokasgrāmata	LV
56	Vartotojo vadovas	LT
59	Kasutusjuhend	ET
61	Посібник користувача	UA
64	Ръководство за потребителя	BG
67	Упатство за употреба	MK
70	Manual do utilizador	PT

A

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Szanowny Kliencie,

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USER MANUAL

Dear Customer,

Thank you for purchasing the Neno Beauty Luvé - a multi-purpose facial device with microcurrent and skin rejuvenation technology. The device is designed for comprehensive and advanced care - it combines up to six effective functions: electroporation, micro-currents, EMS technology, electrical micro-pulses, LED light therapy and sonic vibration. This user manual contains all the information necessary for correct use. Please familiarise yourself with it and keep it in case you need to use it again.

01. PRECAUTIONS

1. Do not use the device on swollen or inflamed skin, or on damaged skin, burned skin or immediately after recent surgery.
2. The device is not recommended for people with electronic implants (e.g. pacemaker, insulin pump).
3. Do not use the device if you have metal components in your mouth (e.g. implants, crowns, braces).
4. Do not use the device if you have active skin diseases (ulcers, psoriasis, eczema, rashes, recurrent herpes), epidermal damage (scars, wounds), hypersensitivity to electromagnetism, pregnancy and lactation, haemophilia/haemorrhagic disorders, epilepsy, atopic dermatitis (AD) and cancer.
5. Do not use the device while taking photosensitising medications (e.g. certain antibiotics, retinoids).
6. Do not use the device in combination with retinol and strong acids - this can lead to excessive skin irritation.
7. The device should not be used by children. Adults with physical, sensory or cognitive limitations should only use the device under the supervision of a carer.

8. If there is any doubt about the suitability of the device for a particular medical condition (e.g. skin diseases, any chronic illnesses, history of cosmetic surgery), consult a doctor.
9. Do not use for too long on one area of skin.
10. Avoid looking directly at the LED light in operation.
11. It is recommended to remove any metal jewellery, such as earrings or chains, before using the device.
12. Avoid using the massage head directly under the eyes, on the eyelids and on the neck in the thyroid region.
13. Always use the lowest intensity level in the eye, nose and mouth area (regardless of the mode selected).
14. Do not use the device longer or more often than recommended. Excessive use can lead to skin irritation.
15. Do not use the device near water or in conditions of high humidity, such as a bath, shower or sauna.
16. If you notice any adverse skin reactions during use, such as redness, itching, burning or other abnormal reactions - discontinue use immediately.
17. The temporary perception of flashing lights during use is a natural phenomenon resulting from the effect of the pulses. However, if it persists afterwards, consult your doctor.
18. A slight tingling sensation felt during use is a normal effect of the pulses. If discomfort occurs, it is recommended to reduce the intensity level.
19. Gentle vibrations or involuntary muscle contractions may occur during use - this is a natural effect of pulse stimulation (microcurrents, EMS, electroporation) and is not a cause for concern.
20. Do not use the device while it is charging. Disconnect the device from the power supply when it is fully charged.
21. Keep the power source away from sources of heat and water while charging.
22. Do not use the product if any component is damaged or missing.
23. Do not disassemble or repair the product yourself.
24. Do not use accessories or spare parts other than those supplied.
25. This product is not a toy. Do not allow children to play with the appliance.
26. Do not store the appliance and its packaging in a place that is easily accessible to children.
27. Do not store the appliance near a heat source and do not expose it directly to sunlight. Do not immerse the device in water.

02. CONTENTS OF THE KIT

1. Massager
2. Storage bag
3. Microfibre cloth
4. Type-C cable
5. Instruction manual

03. DESCRIPTION OF THE DEVICE

SEE FIGURE A

1. Massage head cover
2. Massage head
3. LED light
4. Display
5. Button for changing the intensity of the mode
6. On/off/mode change button
7. USB charging port
8. Type-C cable

04. AIR SHOT

Mode with blue LED light (460 nm), electric micro-pulses and short-wave vibration helps to reduce the appearance of pores, improve skin texture and give skin a smooth and fresh appearance.

Application:

- Wash and dry face thoroughly before use. Mode should only be used on dry skin, without any skincare products.
- The mode offers five levels of intensity that can be adjusted to suit individual skin needs. When using for the first time, it is recommended to start with the lowest level and the intensity can be gradually increased after a few weeks of regular use. For more sensitive areas, such as around the eyes, nose or mouth, always use the lowest intensity level.
- Gently apply the edge or tip of the massage head to the skin - do not use the entire surface of the massage head. Use the mode by making slow, smooth movements along the selected area of the skin or by applying the head directly to areas of concern.
- When using the AIR SHOT mode, you may feel a gentle, pinpoint tingling sensation reminiscent of fine needle-like tingles, as well as smell a slight odour reminiscent of burning fine hairs - this is perfectly normal and safe.
- Do not use AIR SHOT mode on days when you are using retinol or strong acids, especially if you have sensitive skin.
- The mode can be used 2-3 times a week (session maximum 5 minutes).

05. BOOSTER

This mode with yellow LED light (590 nm), gentle vibration and electroporation promotes deep absorption of active ingredients, intensely hydrates and brightens the skin, restoring a healthy glow.

Application:

- Thoroughly wash and dry your face before use (if using after AIR SHOT mode, this step can be skipped), then apply skincare product directly to skin. The best products are those with a gel or water-gel consistency, containing active ingredients such as hyaluronic acid, peptides, niacinamide or ceramides to further help hydrate and regenerate the skin.
- The mode offers five levels of intensity, which can be adjusted to suit individual skin needs. When using for the first time, it is recommended to start with the lowest level, and the intensity can be gradually increased after a few weeks of regular use. For more sensitive areas, such as around the eyes, nose or mouth, always use the lowest intensity level.
- Gently apply the entire surface of the massage head to the skin, using circular motions (from the centre of the face outwards and upwards) to allow the product to be absorbed. If, during the 5-minute BOOSTER mode treatment, you feel that all the skincare product has already been

absorbed, you can apply additional product to ensure smooth gliding of the device and increase comfort.

- The mode can be used 1-3 times a day (5-10 minute session).

06. MC

Mode with green LED light (521 nm), vibration and micro-currents stimulates collagen production, improves skin firmness and smooths fine lines and wrinkles, restoring a youthful appearance to the complexion.

Application:

- Thoroughly wash and dry your face before use (if using after BOOSTER mode, this step can be skipped), then apply the skincare product directly to the skin. For optimal results and comfort, we recommend using the device with a conductive gel suitable for EMS technology or with a light, oil-free facial care product that does not impede the smooth gliding of the device over the skin.
- The mode offers five levels of intensity, which can be adjusted to suit individual skin needs. When using for the first time, it is recommended to start with the lowest level and the intensity can be gradually increased after a few weeks of regular use. For more sensitive areas, such as around the eyes, nose or mouth, always use the lowest intensity level.
- Gently apply the entire surface of the massage head to the skin and move the device slowly across the face, using upward and outward movements - following the natural direction of lymphatic drainage.
- The mode can be used 1-3 times a day (session of 5-10 minutes).

07. DERMA SHOT

Mode with red LED light (630 nm), vibration and EMS technology intensively firms facial muscles, contours and improves skin tone - for a firmer and more naturally sculpted face.

Application:

- Before use, wash and dry your face thoroughly (if using after MC mode, this step can be skipped) and then apply a skincare product directly to the skin. For optimal results and comfort, we recommend using the device with a conductive gel suitable for EMS technology or a light, oil-free facial care product that does not impede the smooth gliding of the device over the skin.
- The mode offers five levels of intensity, which can be adjusted to suit individual skin needs. When using for the first time, it is recommended to start with the lowest level and the intensity can be gradually increased after a few weeks of regular use. For more sensitive areas, such as around the eyes, nose or mouth, always use the lowest intensity level.
- Gently apply the entire surface of the massage head to the skin and move the device slowly across the face in an upward and outward motion - following the natural direction of lymphatic drainage. In selected areas, you can pause the device for 1-2 seconds to concentrate the pulses on muscle tension points, which will help to relax them.
- Before using the DERMA SHOT mode, it is recommended to use the MC mode first.
- The mode can be used 1-3 times a day (session of 5-10 minutes).

08. OPERATION OF THE DEVICE

1. To switch on the device, press and hold button A6 (a beep will sound). Use the same button to set one of the four available modes and then use button A5 to set one of the five intensity levels (the device will emit a single beep with each change).
2. The device has a built-in 5-minute timer and a skin contact sensor. The countdown timer only starts when the massage head touches the skin - otherwise the timer stops. In addition, the LED light only switches on when it comes into contact with the skin, which increases user safety as it eliminates the risk of accidental light emission into the air or eyes.
3. After 5 minutes, the device automatically goes into a resting state for 9 seconds. During this time, the A6 button can be pressed to select the next mode of operation. If no operation is performed during these 9 seconds, the device switches off automatically.
4. If the operating mode is changed during the treatment, the countdown starts from the beginning. In contrast, changing the intensity in the same mode does not reset the timer - the time adds up so that the user still receives the full 5 minutes of treatment.
5. The device allows all four modes to be used in a single treatment cycle. However, they must be used in the recommended order for best results and safety: AIR SHOT -> BOOSTER -> MC -> DERMA SHOT. If you use several modes in the correct order, there is no need to wash and dry your face every time you change modes - just apply the appropriate skincare product or conductive gel. Remember that only the AIR SHOT mode should be used on dry skin, without using a skincare product.
6. The device operates in pulsed mode - a deliberate and safe action that allows for effective skin stimulation.
7. To switch off the device before 5 minutes have elapsed, press and hold button A6 - a double beep will sound.
8. Once the massage has finished, you can move on to the next stages of care.

09. CLEANING AND STORAGE

1. After each use, wipe the massage head with a dry or slightly damp cloth, taking care to remove any residual cosmetics used during the massage.
2. If necessary, wipe the head dry with a clean cloth and do not use harsh detergents or alcohol.
3. Do not immerse the device in water and do not use highly corrosive detergents or alcohol.
4. When the product is not in use, store it in the bag provided, in a dry, cool place, away from direct sunlight and heat sources.

10. CHARGING THE DEVICE

1. When charging, the charging icon is visible on the screen and the green ring shows the charging progress.
2. When the device is fully charged, the charging ring lights up fully green.
3. When the battery level is low, the battery level indicator on the display turns red and displays one dash.
4. When the device is discharged, the battery level indicator flashes three times, the device beeps 3 times and automatically switches off.

11. SPECIFICATIONS

Power supply: 5V 0.7A

Power: 5W

Battery: Lithium-polymer 2200mAh

Timer: 5min

Material: ABS + PP + 304 stainless steel

Dimensions: 38.2x182.7x 38.6 mm
Weight: 304g

12. WARRANTY CARD

The product is covered by a 24-month warranty. Warranty conditions can be found at: <https://neno.pl/gwarancja>
The details, contact, and service address can be found at: <https://neno.pl/kontakt>
Specifications and contents are subject to change without notice. We apologise for any inconvenience.

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