nenob[®] beauty

Naya


Table of contents

Spis treści

4	Instrukcja obsługi	PL
5	User manual	EN
7	Bedienungsanleitung	DE
9	Manuel de l'utilisateur	FR
10	Návod k použití	CZ
12	Používateľská príručka	SK
14	Használati utasítás	HU
15	Användarhandbok	SE
17	Käyttäjän käsikirja	FI
18	Brakerhåndbok	NO
20	Brugermanual	DK
22	Gebruikershandleiding	NL
23	Manual del usuario	ES
25	Manuale utente	IT
27	Manual de utilizare	RO
28	Korisnički priručnik	HR/BA
30	Uputstvo za upotrebu	RS/ME
31	Uporabniški pravilnik	SI
33	Εγχειρίδιο χρήστη	GR
35	Lietotāja rokasgrāmata	LV
36	Vartotojo vadovas	LT
38	Kasutusjuhend	ET
40	Посібник користувача	UA
41	Ръководство за потребителя	BG
43	Упатство за употреба	MK
45	Manual do utilizador	PT

A


1

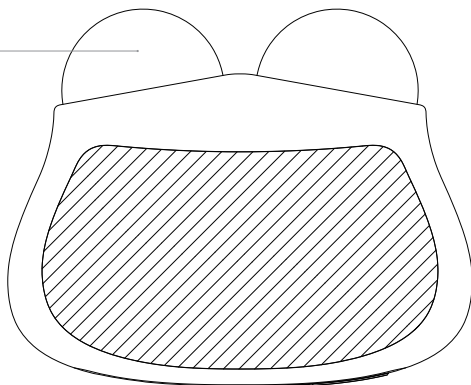
4

3

2

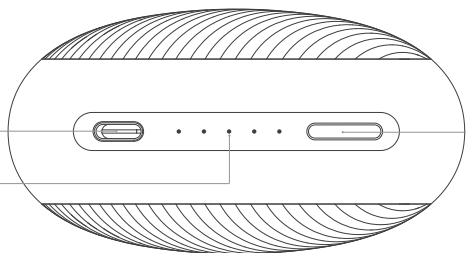5

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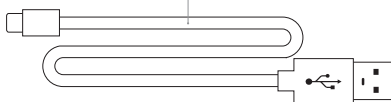Szanowny Kliencie,

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USER MANUAL

Dear Customer,

Thank you for purchasing the Neno Beauty Naya facial massager with EMS technology and thermotherapy function. The device is designed for daily facial care - it uses three effective functions to promote skin firming, improve microcirculation and help contour the face, restoring a healthy and radiant appearance. This user manual contains all the information necessary for correct use. Please familiarise yourself with it and keep it in case you need to use it again.

01. PRECAUTIONS

1. Do not use the device on swollen or inflamed skin, or on damaged skin, burned skin or immediately after recent surgery.
2. The device is not recommended for people with electronic implants (e.g. pacemaker, insulin pump).
3. Do not use the device if you have metal components in your mouth (e.g. implants, crowns, braces).
4. The device should not be used in cases of active skin diseases (ulcers, psoriasis, eczema, rashes, recurrent herpes), epidermal damage (scars, wounds), hypersensitivity to electromagnetic fields, pregnancy and lactation, haemophilia/haemorrhagic disorders, epilepsy, atopic dermatitis (AD) and cancer.
5. Do not use the device while taking photosensitising drugs (e.g. certain antibiotics, retinoids).

6. The device should not be used by children. Adults with physical, sensory or cognitive limitations should only use the device under adult supervision.
7. If you have any doubts about the suitability of the device for a specific medical condition (e.g. skin diseases, any chronic illnesses, history of cosmetic surgery), consult your doctor.
8. Do not apply too long on one area of skin.
9. It is recommended to remove any metal jewellery, such as earrings or chains, before using the device.
10. Avoid using massage heads around the eyes and on the neck in the thyroid region.
11. Do not use the device longer or more often than recommended. Excessive use can lead to skin irritation.
12. Do not use the device near water or in high humidity conditions such as a bath, shower or sauna.
13. If you notice any adverse skin reactions during use, such as redness, itching, burning or other abnormal reactions - discontinue use immediately.
14. The temporary perception of flashing lights during use is a natural phenomenon resulting from the effects of the pulses. However, if it persists afterwards, consult your doctor.
15. The slight tingling sensation felt when using the device is a normal effect of the EMS pulses. In case of discomfort, it is recommended to reduce the intensity level.
16. Do not use the device while it is charging. Disconnect the device from the power supply when it is fully charged.
17. Keep the power source away from sources of heat and water while charging.
18. Do not use the product if any component is damaged or missing.
19. Do not disassemble or repair the product yourself.
20. Do not use accessories or spare parts other than those supplied.
21. This product is not a toy. Do not allow children to play with the device.
22. Do not store the device and its packaging in a place that is easily accessible to children.
23. Do not store the device near a heat source and do not expose it directly to sunlight. Do not immerse the device in water.

02. CONTENTS OF THE KIT

1. Massager
2. Base
3. Microfibre cloth
4. Type-C cable
5. Operating instructions

03. DESCRIPTION OF THE DEVICE

SEE FIGURE A

1. Double massage head
2. Button for switching on/off/changing the intensity level
3. Intensity level indicator / Power indicator
4. USB charging port
5. Type-C cable

04. OPERATION OF THE DEVICE

1. The device combines the action of three functions: EMS (Electrical Muscle Stimulation), thermotherapy and vibration massage, offering 5 levels of intensity that vary in EMS stimulation range and vibration strength. This makes it easy to adapt the action of the device to the user's individual preferences and skin condition.
2. Wash and dry your face thoroughly before use. For optimal results and comfort, we recommend using the device with a conductive gel suitable for EMS technology or with a light, oil-free facial cosmetic that does not impede the smooth gliding of the device over the skin.
3. To switch the unit on, press and hold the A2 button (a beep will sound). Then use the same button to set one of the five intensity levels - the unit will beep once each time the level changes. The intensity level indicator lights will illuminate in the number corresponding to the selected level (1-5). The device works in pulsed mode - a targeted and safe action that allows for effective skin stimulation.
4. During the massage, move the device slowly across the face, making upward and outward movements - following the natural direction of lymphatic drainage. The device signals the change of area every 15 seconds (beeps twice).
5. The device has a built-in 6-minute timer: the device will switch off six minutes after one of the levels has been set. The time does not add up when changing intensity levels during their operation.
6. To switch off the unit before 6 minutes have elapsed, press and hold the A2 button - a double beep will sound.
7. Once the massage is over, you can move on to the next stages of care.
8. The device has a memory mode - when the massager is switched on again, it will automatically start in the last mode used.
9. Use the massager 2-4 times a week (for 6 minutes), choosing an intensity level that suits your needs.

05. CLEANING AND STORAGE

1. Wipe the head with a dry or slightly damp cloth after each use.
2. Leave to dry in a dry place, away from direct sunlight and heat sources.
3. Do not immerse the device in water and do not use highly corrosive detergents or alcohol.
4. When not in use, the product should be stored in a dry place on the supplied stand.

06. CHARGING THE DEVICE

1. During charging, the power indicator LEDs flash successively from left to right. The device cannot be used during charging.
2. The power indicator lights illuminate continuously for 5 minutes when the unit is fully charged and then go out.
3. When the unit is switched off, press the A2 button briefly to check the battery level - this will be shown by the number of LEDs lit. One LED indicates approximately 25% charge level. The LEDs will go out after 2 seconds.
4. When the battery level is low, the power indicator LEDs flash and the unit will beep twice to indicate that the unit needs to be recharged.

5. If an attempt is made to start the unit with a fully discharged battery, a beep will sound three times, the power indicator LEDs will flash three times and the unit will automatically switch off.

07. SPECIFICATION

Power supply: 5V 0.7A

Power: 6W

Battery: Lithium-polymer 1000mAh

Timer: 6min

Material: ABS + PC + silicone + zinc alloy

Dimensions: 84x45x 68 mm

Weight: 172g

08. WARRANTY CARD

The product comes with a 24-month warranty. Warranty terms and conditions can be found at: <https://neno.pl/gwarancja>

Details, contact and service address can be found at: <https://neno.pl/kontakt>

Specifications and contents are subject to change without notice. We apologise for any inconvenience.

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