



MECHANICZNA WAGA OSOBOWA
MECHANICAL PERSONAL SCALE

HELIA

Webber⁺
med

PL

INSTRUKCJA OBSŁUGI

model: **MBS-605T**

Gratulacje, zakupiłeś wagę mechaniczną Webbermed wykonaną z materiałów wysokiej jakości.

Przed jej użyciem, prosimy o zapoznanie się z instrukcją obsługi.

INSTRUKCJA UŻYTKOWANIA



Wskaźnik BMI

Wskaźnik masy

Skala masy

Obrotowe ustawienie wzrostu użytkownika

1. Umieść wagę na płaskiej powierzchni.
2. Przed użyciem ustaw wskaźnik masy na zero.
3. Obróć pierścień, aby ustawić swój wzrost na poziomie zera.
4. Wejdź na wagę, aby odczytać masę i wskaźnik BMI.

PIELĘGNACJA URZĄDZENIA


- Stań delikatnie, wskaźnik może się zmieniać, jeśli się poruszyysz.
- Nie używaj ściernych środków czyszczących.
- Nie używaj wagi do celów komercyjnych.



Congratulations, you have purchased one of the Webbermed mechanical personal scale made of high quality materials.

Please read usage description before use it.

OPERATION GUIDANCE

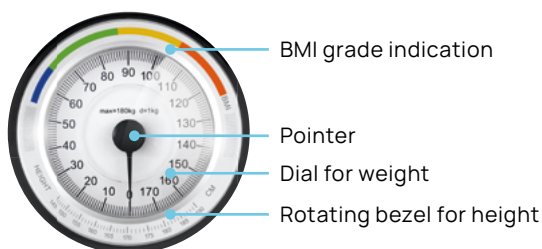

1. Place your scale on flat surface.
2. Adjust dial to zero before use.
3. Rotate the bezel to set your height aligned with zero.
4. Step on scale, the scale indicates the weight & BMI grade.

CARE OF YOUR SCALE

- Stand gently, display will vary if you move.
- Do not use abrasive cleaners.
- Do not use for commercial purpose.

