BEAUTIFLY

Aura NovaSkin

PRO

series

LED Mask
User manual

/beautifuly.global

/beautifuly_global


/@beautifuly.global


/@beautifuly.global


f1

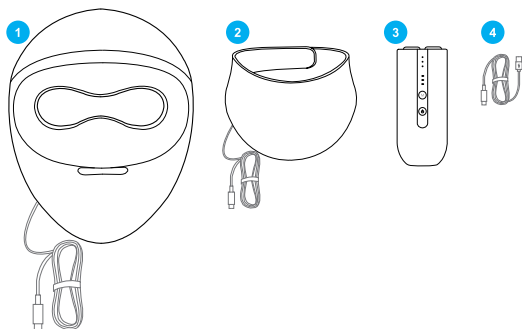f3

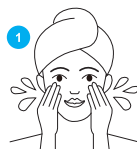f2

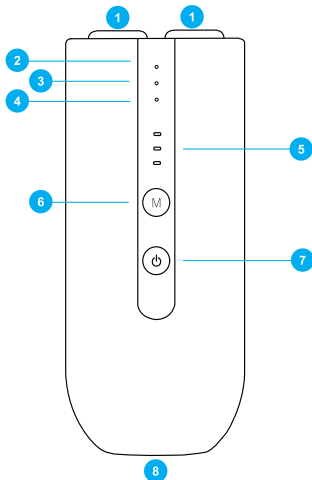
f4

facebook.com/beautifly.global

youtube.com/@beautifly.global

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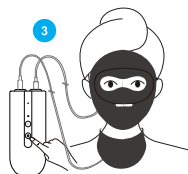Operating Voltage:	3.7V
Charging Voltage:	DC5V 1A
Number of LEDs:	Face – 444 pcs / Neck – 58 pcs
Battery Capacity:	2000 mAh
Material:	ABS Silicone
Red Light + NIR:	630 ± 10nm + 850 ± 5 nm
Green Light:	505 ± 5nm
Blue Light:	470nm

Safety:

Read this user manual before using this product and follow the instructions contained therein when using the device.

Kit Contents [page 2-f1]:

1. Face Mask / 2. Neck Mask / 3. Remote Control / 4. USB Cable

Remote Control Description [page 2-f2]:

1. Power Output x2 / 2. IR + Red Light Mode / 3. Green Light Mode / 4. Blue Light mode /
5. Mode Switch Button / 6. Intensity Level (Low, Medium, High) / 7. ON/OFF Button / Intensity Adjustment / 8. Charging Port.

Device Operation:

1. Turning On/Off: Hold the button for 3 seconds to turn the device on or off.
2. Automatic Mode: Once turned on, the device enters automatic mode and begins a 20-minute phototherapy cycle consisting of three stages:
- Stage 1 (10 min) – Infrared (invisible) light + red light; improved circulation and collagen stimulation.
- Stage 2 (5 min) – Green light; acne reduction and skin cleansing.
- Stage 3 (5 min) – Blue light; antibacterial properties, helps alleviate acne discoloration.
The automated program allows for a convenient, complete phototherapy session in the comfort of your home.
3. Adjusting the intensity: During operation, briefly press the "ON/OFF" button to change the brightness (L / M / H).
4. Charging: A flashing battery indicator indicates low battery. Plug the device into a power source as soon as possible to avoid interrupting your treatment.
5. Recommended use: Use the rejuvenating mask 3–5 times a week.

How to use the device [page 2-f2]:

1. Remove makeup and thoroughly cleanse your face.
2. Apply a mask or serum, then secure the device to your head using the Velcro straps.
3. Hold the power button for 3 seconds to turn on the device. The program runs automatically for 20 minutes.
4. Remove the mask and apply moisturizer.

Device Care and Maintenance:

1. Store this product in a dry place at room temperature and avoid exposure to dust, direct sunlight, or humid conditions.
2. Clean the mask (mask only, excluding electronic parts) before and after use. Note: The controller and cables should not come into contact with water. After use, ensure you unplug the controller from the mask and from the power source.
3. This product contains a lithium polymer battery. When disposing of the device, please take it to a local battery recycling point. It is recommended to separate it from other household waste to protect the environment.

Precautions:

Keep the product out of the reach of small children. Handle the product with caution during use. It should not be used by pregnant women, those planning to become pregnant, breastfeeding, or people sensitive to light. Do not use it while charging. Do not apply excessive pressure to the device to avoid deformation and malfunction. Do not repair, disassemble, or modify the device yourself. Do not drop or hit the device to prevent malfunction. Note: The remote control and cables should not come into contact with water. After each treatment, wipe the product with a damp cloth. Do not wash the device under water or submerge it. Do not use abrasive substances to clean the device. Store it in a cool, dry place. The device has an automatic shut-off function after 20 minutes of use. Regardless of the safety precautions taken, the device may cause skin redness or irritation. In such cases, discontinue use immediately and consult a doctor. If the device is used more frequently than recommended, there is a risk of skin irritation or inflammation. Do not use a solarium or expose the skin to sunlight for at least 48 hours before using the device and for at least 48 hours after using it. The device is not suitable for use on individuals under 18 years of age. Do not use on areas of skin where there is an increased risk of skin cancer, on open wounds or rashes, on swollen, red, irritated, infected, or inflamed areas of skin, or on skin eruptions. Wear eye protection when using IPL or laser devices to avoid exposure to direct or reflected radiation. Before using the device, it is recommended to consult a physician and, if necessary, perform diagnostic testing of the area being treated. Treatment should not be performed until documentation of such consultation is available.

More information about the company and products [page 2-f4]:

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[EN] Rules for handling electrical and electronic equipment: Used electrical and electronic equipment, batteries, and accumulators should not be disposed of with other waste. This is indicated by the crossed-out waste bin symbol on the product. Electrical and electronic equipment, batteries, and accumulators may contain accessories, components, and parts that are harmful to the environment and human health. Improper handling or their consequences may result in damage to the environment or the environment. Proper household handling of used electrical and electronic equipment, batteries, and accumulators is essential for the proper disposal of devices and equipment, including used electrical and electronic equipment. Used electrical and electronic equipment, batteries, and accumulators can be disposed of separately. The distributor is obligated to collect household equipment free of charge at the point of sale, provided the used equipment is of the same type and performs the same functions as commercial equipment. A distributor delivering household appliances to a purchaser is obligated to collect waste household appliances free of charge at the point of delivery, provided the waste appliances are of the same type and have served the same function as the appliances. A distributor operating a retail business with an online sales area of 400 m² of household appliances designated for sale pursuant to Article 2, Section 19 of the Act of 27 March 2003 on Spatial Planning and Development is entitled to free membership in the same entity or its immediate vicinity for discarded household appliances with external dimensions not exceeding 25 cm, without the need to purchase new household appliances. The product was registered on the market after August 13, 2005.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product Aura NovaSkin Pro complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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The full text of the declaration of conformity is available at the following internet address: www.beautify.eu/deklaracje
Importer: Beautify Sp. z o.o. Plac Bankowy 2, 00-095 Warszawa, Polska

Model: Aura NovaSkin Pro
SN: 05265903003703297
Marka: Beautify
e-mail: kontakt@beautify.eu