BEAUTIFLY

SkinWOW Shot
Golden
series

LED light
therapy device

User manual

/beautifuly.global

/beautifuly_global

@beautifuly.global

@beautifuly.global

f1


facebook.com/beautifly.global

youtube.com/@beautifly.global

f2

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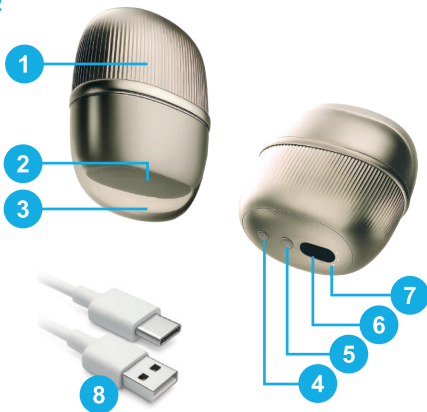Device description (see page 2, f2)

1. Main unit
2. Matrix probe
3. Dust cover
4. Power / light switch button
5. Intensity adjustment button
6. Charging port (USB-C)
7. Intensity indicator / charging indicator
8. USB cable

Technical specifications:

- Charging type: USB-C
 Charging time: 3 hours
 Input voltage: 5V
 Battery capacity: 1100 mAh
 Product power: ≤10W
 Frequency: 1000KHz (1MHz)
 Number of intensity levels: 3
 Number of light color levels: 3
 Set contents: 1x main unit, 1x dust cover, 1x USB cable,
 1x user manual, 1x warranty card

The device has three modes:

1. Stamping Mode (red light 650 nm)—firms the skin, supports collagen regeneration, prevents wrinkle formation
2. EMS Mode (blue light 465 nm)—reduces skin puffiness, lifts
3. Gold Slide Lift Mode (orange light 600 nm)—reduces wrinkles, improves skin laxity, increases dermal thickness and density

Device operation

1. To turn on the device, hold the power button. Then you can select the desired working mode by toggling through the modes.
2. After switching on, you can adjust the intensity using the intensity adjustment button.
3. To switch off, hold the power button. Protective auto-off function: If no button operation is performed for 5 minutes, the device will switch off automatically after 5 minutes of continuous use.
4. To turn on vibration, double-click the power button. Another double-click switches vibration off. Changing the mode also automatically turns vibration off. Regardless of the selected mode, vibration is disabled by default.

Preparation for use

1. Remove make-up and thoroughly cleanse face and neck to ensure the skin is clean.
2. You can apply gel to the face or squeeze it directly onto the device probe.
3. A water-based cosmetic increases RF conductivity and improves penetration of active ingredients.
3. Switch on the device and select the desired mode for individual use, or use each of the 3 modes once.
4. After the session, wash the face with water and apply an appropriate amount of skincare products to moisturise as needed.
5. After the session use a slightly damp, soft cloth to wipe residues from the beauty device. If there is a lot of residue, you can rinse the device directly with water, then dry it with a dry cloth.

Recommended usage frequency

1. Stamping Mode (red light): Normal skin—3 times a week for 10-15 min, Dry skin—2 times a week for 10-15 min, Oily skin—3-4 times a week for 15-20 min, Sensitive skin—once a week for 5-10 min.
2. EMS Mode (blue light): Use every 2 days for 5 min.
3. Gold Slide Lift Mode (orange light): Use every 2 days for 5 min.

Charging

1. Before charging, make sure the device's charging port is dry and free of water and skincare residues.
2. Connect the USB charging cable to a 5V USB charger (you may also use a computer USB port or power bank), then insert the Type-C plug into the port in the charging base.
3. During charging the indicator will pulse. When fully charged, the indicator will stay steadily lit, no longer pulsing or flashing.
4. Do not switch the device on while it is connected to the charger.
5. Fully charge the device at least once every two months to protect battery life.

Cleaning and maintenance

1. Ensure the device is switched off before cleaning.
2. Do not use aggressive cleaning agents, petrol, acetone, etc. to clean the housing.
3. To maintain treatment effectiveness, clean the probe after each use.
4. Regular cleaning and proper maintenance ensure optimal performance and extend the device's lifespan.

1. Before using the purchased product, read the product manual and these safety warnings, maintenance and storage instructions, as well as rules for handling waste electrical and electronic equipment. Follow the recommendations contained therein and use the product as intended.

2. Users and consumers are advised to consult a physician and undergo any necessary diagnostic tests of the body parts to be treated before using the device. The manufacturer recommends not undergoing a session before obtaining documentation from such a consultation.
3. Excessive use (more than 3 sessions per day on the same spot) may cause irritation or thermal burns.
4. Maintain at least a 4-hour break between sessions on the same skin area.
5. Persons under 18 years of age must not use the device.
6. Do not use around the eyeball, mouth, throat.
7. Do not use on skin areas with an increased likelihood of skin cancer, on open wounds, scars or rashes, on swollen, reddened, irritated, infected or inflamed skin areas.
8. Persons suffering from blood diseases, tuberculosis, epilepsy, diabetes, coeliac disease, malignant tumours, hypertension, fever must not use the device.
9. Do not use if you are pregnant, menstruating or breastfeeding.
10. The device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or persons lacking experience and knowledge, unless they are supervised or instructed by a person responsible for their safety.
11. It is not recommended to use the device after prolonged skin exposure to sunlight or tanning beds.
12. Eye protection is recommended for persons exposed to radiation due to reflection, misuse or unskilled handling of the device during a session.
13. Do not use too long on the same skin area.
14. The device must not be used if you have electronic devices inside the body, e.g. a pacemaker.
15. Remove watches, glasses, jewellery and other metal objects before using the device.
16. Do not use the product if it does not work or works incorrectly, or has been dropped or otherwise damaged. In such a situation, do not open or repair the product yourself. If the product is still under statutory warranty or guarantee, contact the seller or guarantor. If the warranty period has expired, or the defect is not inherent to the product, take it to an appropriate service centre.
17. Do not use the product near sources of heat, water or moisture, and do not immerse it in water or other liquids.
18. The product is intended for use at room temperature, in the range of 10-40 °C.
19. Do not expose the product to direct sunlight.
20. After use, store the product in a dry place at room temperature, in the range of 10-40 °C, out of reach of children and animals.
21. If the product becomes dirty, disconnect it from the power source, wipe with a damp cloth and allow it to dry before storage. Do not use abrasive materials or caustic cleaning agents, as this may damage the product.
22. Do not use the product if you notice a damaged power cord or plug, an unpleasant smell (e.g. burning), overheating or smoke. In such a situation, do not open or repair the product yourself. If the product is still under statutory warranty or guarantee contact the seller or guarantor. If the warranty period has expired, or the defect is not inherent to the product, take it to an appropriate service centre.
23. Do not leave the product unattended while it is switched on.
24. The power source must match the product's required voltage.
25. After charging, disconnect the product from the power source and allow it to cool before moving or storing it.
26. Do not twist or coil the power cable, as this may damage it.
27. Accessories for the product that do not originate from the manufacturer may lead to malfunction or damage.
28. Consult a dermatologist before using the device if you have any skin health problems.
29. If after using the device you experience skin discomfort (tightness, discomfort, irritation, redness, burning, etc.), notice unexpected side effects or prolonged redness, stop using immediately and consult a physician. Temporary, mild redness is normal, however.
30. Do not look directly at the probe's light source for an extended time, as this may damage your eyesight.
31. After the declared service life (2 years) has elapsed, take the device to an e-waste collection point.
32. Do not share the device with other people, as this may lead to cross-infection by bacteria.

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[EN] Used electrical and electronic equipment, as well as batteries and accumulators, must not be mixed with other waste. This is indicated by the symbol of a crossed-out bin placed on the product packaging. Electrical and electronic equipment, as well as batteries and accumulators, may contain hazardous substances, mixtures, and components harmful to the environment and human health. Improper handling or damage to them can lead to health risks or environmental pollution. Proper handling of used electrical and electronic equipment, as well as batteries and accumulators, in households positively contributes to reuse and recovery, including recycling of the used equipment. Used electrical and electronic devices, as well as batteries and accumulators, can be returned at local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and served the same functions as the equipment sold. When delivering equipment intended for households, the distributor is obliged to collect used equipment from households free of charge at the delivery location, provided that the used equipment is of the same type and served the same functions as the delivered equipment. A distributor operating a retail unit with a sales area, as defined in Article 2, point 19 of the Act of March 27, 2003, on spatial planning and development, of at least 400 m² dedicated to the sale of equipment intended for households, is obliged to accept used equipment from households free of charge at this unit or in its immediate vicinity, provided that none of the external dimensions exceeds 25 cm, without the necessity of purchasing new equipment intended for households. The product was placed on the market after August 13, 2005.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product SkinWOW Shot Golden complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances. The full text of the declaration of conformity is available at the following internet address:

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