

FOOT MASSAGER MANUAL

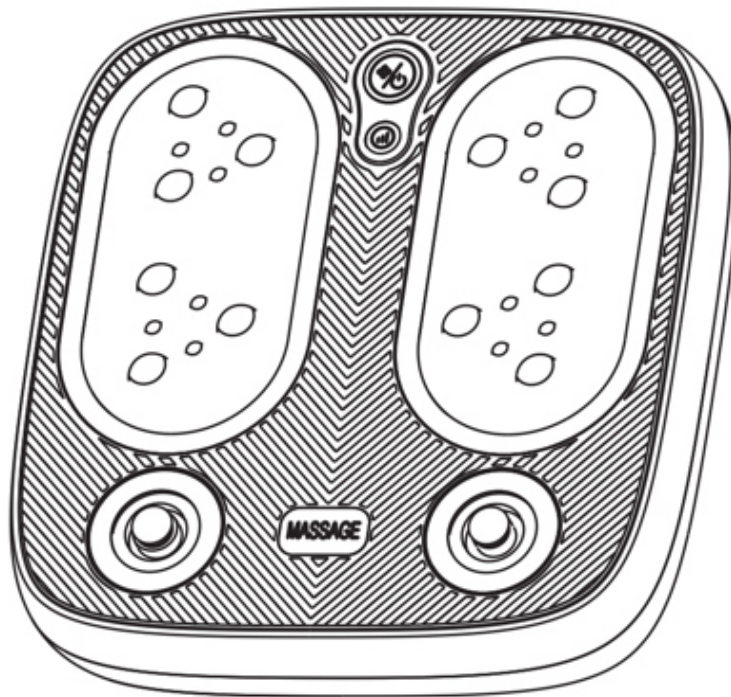DIFORO
DiF-Feet

Thank you for purchasing this air pressure massager, please read this manual thoroughly before using it, please pay attention to the safety precaution, keep this manual in a safe place for later review.

1.Overview

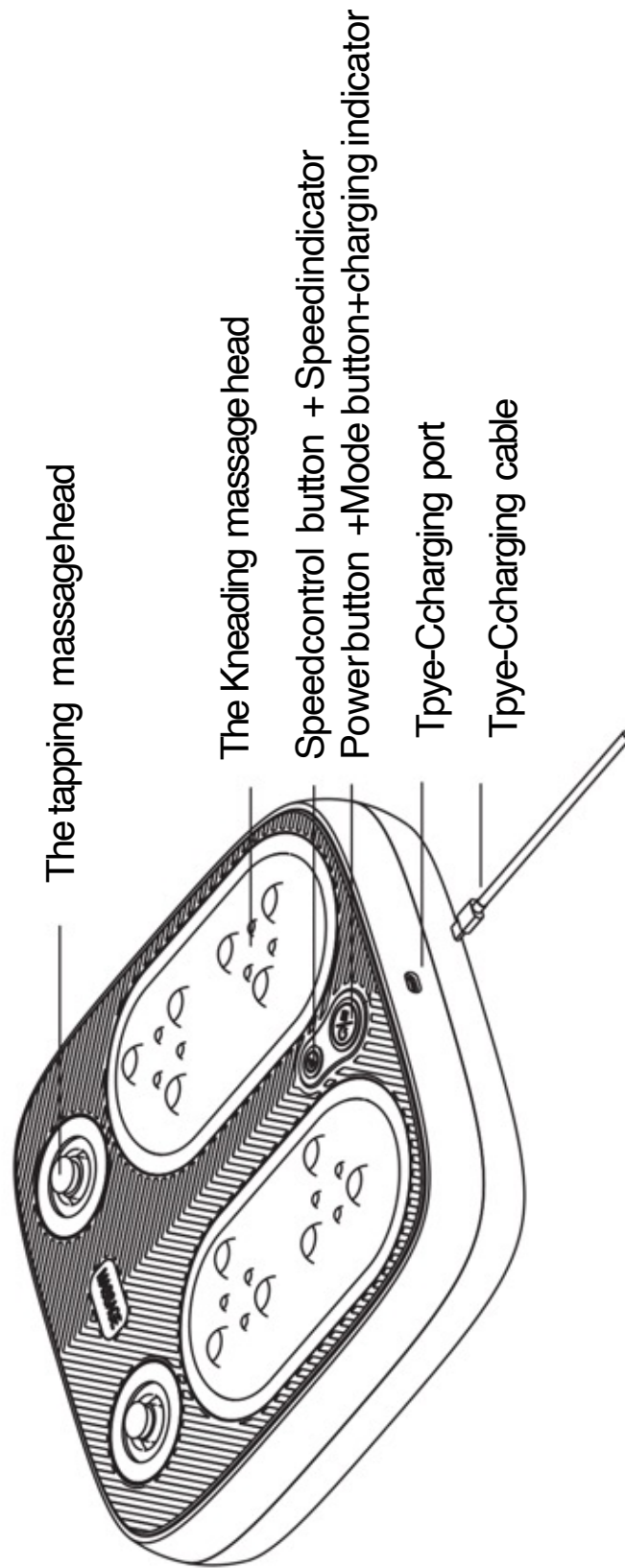Do not use this product while standing

List: DiF-Feet ,Type-C charging cable,manual

2. Operating instructions

- 1. Power/Mode button: Long press to activate, the short press switch mode, when turn on the running indicator light that stays on (blue light), default mode 1, press once to switch mode 2, then press again to switch mode 3 in this loop.

Three Modes: Kneading + tapping mode (without heating); Kneading + tapping mode (with heating); Kneading + tapping mode (intermittent rhythm without heating).

- 2. Speed button: After opening, press to adjust the kneading and tapping speed, Enable default speed mode 1 (indicator light is blue); Press once to switch to speed 2 (the indicator light is orange); Press again to switch to speed 3 (the indicator light is red), and repeat the cycle of switching. Massage head rotation speed: L1-24 rpm, L2-27 rpm, L3-28 rpm; Tapping (motor) rotation speed: L1-2340 rpm, L2-2520 rpm, L3-2700 rpm.

- 3. Heating: Switches to mode 2, the red light heating is working. In other modes, heating is not activated.

3. Blockage protection

After 5 seconds of blockage, the massage head stops running for 3 seconds and the machine automatically continues to run.

4. Intelligent protection

This machine has intelligent protection function. To prevent damage during prolonged massage, enter the closed state after working continuously for 10 minutes..

To extend battery life, the machine cannot be turned on after being turned off due to low battery. It needs to be charged before the machine can be turned on and used.

When using, do not forcefully press the massage head to obstruct the motor rotation, otherwise it will reduce the service life of the equipment! When the battery is low, continuing the operation will result in automatic shutdown.

5.Charging

When using the Type-C charging cable for charging, the power mode button indicator light that stays on blue breathing light. When fully charged, the power mode button indicator light stays on blue constant light. When the power button indicator light flashes, the product enters a low battery state. Please charge it in a timely manner.

6.Parameters

Product Name: Tapping and Kneading Foot Massager

Model: DiF-Feet

Battery capacity: 2500mAh

Charging input: 5V 2A

Power: 24W

Product weight: 1.92kgs

Product size: 318x304x75mm

Heating temperature: About 45 degrees

Taboo

1. Prohibit cancer patients from using it
2. It is prohibited for patients who have undergone metal bone needle internal fixation surgery to use it directly at the fracture site
3. Prohibit patients with osteoporosis from using it
4. Prohibit use in patients with venous thrombosis and pulsation
5. Prohibited for hypertensive patients to use
6. Pregnant women are prohibited from using it
7. Prohibit patients with venous dilation from using it
8. Prohibited for use by patients with acute inflammation
9. Do not massage the bone area
10. Prohibited from use in postoperative areas
11. Massage the same position for no more than 20 minutes.

Safety instructions

1. Keep the massage area dry and clean. It is recommended to use it with shoe covers or socks, and gently press and move it.
2. High speed and high pressure can only be used on soft tissue parts of the feet without causing pain or discomfort, and cannot be used on any bone parts.
3. Frequent high-speed and high-pressure pressure on the same area can cause abrasions, pain, and discomfort. Stop using immediately.
4. Fingers, hair, or any part of the body hair should be kept away from the stretching part of the massage head to prevent pinching pain and injury.
5. Do not discard or misuse this product.
6. Please do not make any modifications to this product without authorization.
7. Do not operate or charge this product in an unmanned state.
8. This product is only for adult use.
9. For individuals with contraindications and symptoms, please use this product with the permission of a doctor.

10. The appliance has a heating surface. Attention must be paid when using it for people who are not sensitive to heat

Pay attention to

1. Please ensure that there are no obstacles in front of the massage head when starting this product. If it is obstructed by external forces, it will pause briefly and then try starting again.
2. This product is not waterproof, please do not let it come into contact with water or other liquids.
3. Do not disassemble the entire product, If accidental disassembly causes abnormal operation, please send it back for re debugging. All expenses incurred as a result shall be borne by the customer.
4. Long term continuous operation may cause the surface temperature of the motor to rise, please operate with caution.
5. It is strictly prohibited to drop or collide this product, as it may cause damage and result in abnormal operation. If there is any malfunction caused by falling or collision, please contact after-sales service in a timely manner.

Storage and maintenance

1. Please store this product out of reach of children and pets.
2. It is prohibited to place this product near heat sources (such as stoves or heaters) or in cars on hot days.
3. Please store in a dry environment.
4. Do not overcharge or over discharge the battery, otherwise it will cause damage to the battery cells
5. Do not use this product in environments where the temperature is too high or too low.

People unsuitable to use

1. This product has no fireproof function, please keep away from fire or corrosive liquid;
2. Do not disassemble or repair this product by yourself;
3. Please keep it at a place out of reach of children;
4. This product has waterproof function;
5. Do not press it heavily so as to avoid deformation and failure.

Warranty

The manufacturer guarantees its manufactured products for a period of 24 months from the date of purchase. Accessories included with the original product are warranted for a period of 12 months from the date of purchase.

For the efficient operation of products and complaints, we recommend.

- ① Study and follow the operating instructions carefully.
- ② Keep the manual, warranty card, proof of purchase.
- ③ Use the original accessories supplied with the device.
- ④ In case of problems with the device, contact the Manufacturer.
- ⑤ Do not make repairs or modifications yourself.
- ⑥ Use the product according to its intended use.

During the warranty period, all repairs are performed free of charge.

Any defects or faults will be removed within a maximum period of 14 days, calculated from the time of delivery of the device to the manufacturer's authorized service center.

Limitations

The warranty will not be recognized, in case of mechanical damage due to inattention or fortuitous events.

If the device is used contrary to its intended use and specifications.

Batteries are consumables, the warranty covers 6 months. Batteries cannot be replaced outside of an authorized service center.

Contact details

DIFORO Brand

CUSTOMER SERVICE
+48 453 271 100

Consultation hours : 08:00 ~ 16:00
(Monday to Friday)

PRODUCENT

DIFORO by VOBÉ sp. z o.o.
Rudna Mała 167A, 36-054 Mrowla, POLAND
www.diforo.pl get@diforo.pl
NIP: 5252799270, KRS 0000800189


Warranty card

Product Name	Beauty LED Mask	Model	
Distributor information			
Buyer information	Name:		
	Address:		
	Tel:		
Cause of failure			

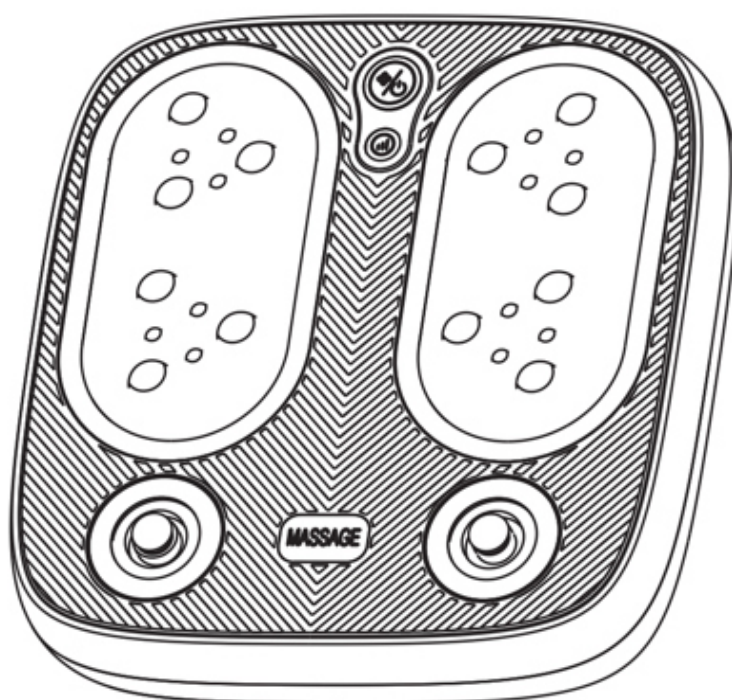Please fill out the warranty card in detail when repairing. Thank you for your cooperation!

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