
EN

DIFORO

Beauty LED Mask

DiF-Eluminna

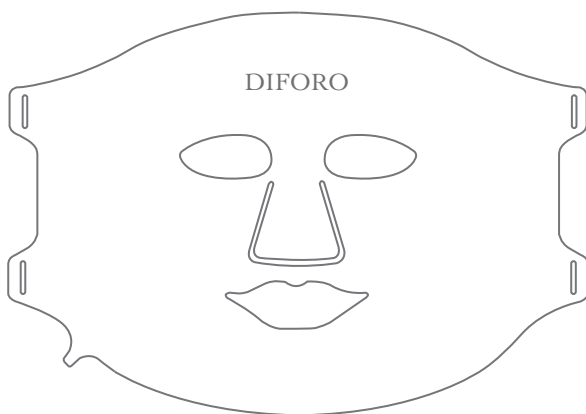Vobe sp. z o.o., Rudna Mała 167A

36-054 Mrowla, Poland

Tel: +48 453 271 100

e-mail: get@diforo.pl

Please read this entire manual before using

Thank you for purchasing the spectrum skin rejuvenation mask, please read this manual carefully before using the mask. Please pay special attention to the safety matters so as to operate and use the mask correctly.

This manual provides users with safe and effective ways of operation. After reading, please keep it properly for future reference.


Content

- Technology Introduction/product list -----
- Controller -----
- Quiiuck Start Guide -----
- Operation / people unsuitable to use -----
- Product specification/troubleshooting -----
- Warranty card -----


Technology Introduction

Wavelength of product.
NIR: 850nm
red light: 630nm
blue light: 415nm
yellow light: 605nm



Product list

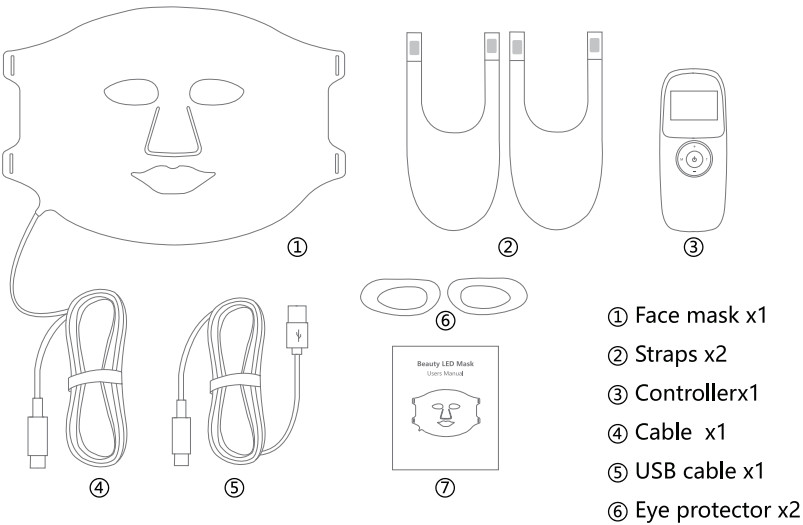Controller

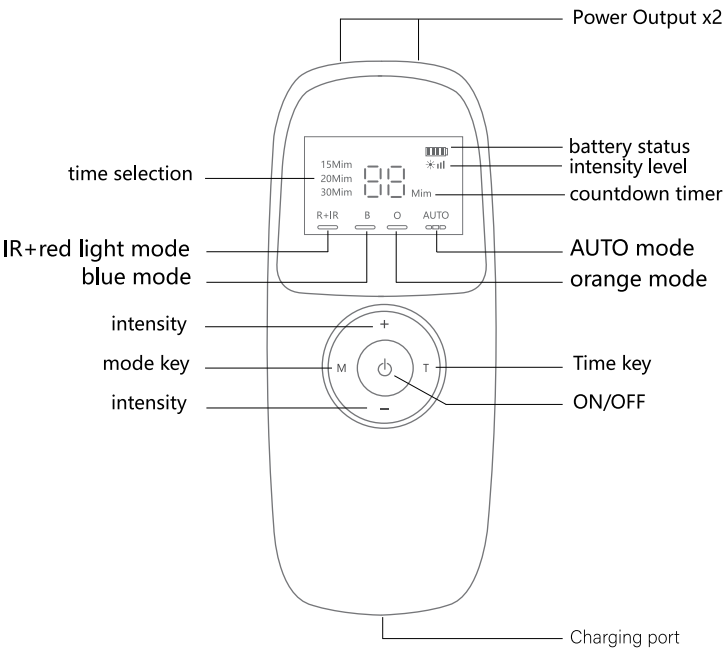Quiick Start Guide

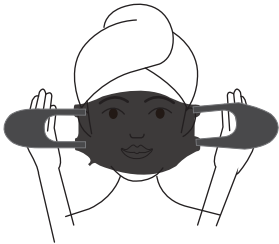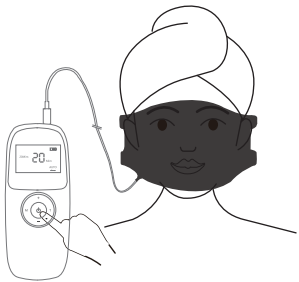 step 1: Remove makeup and clean your face.	 step 2: Use with a mask or serum Secure the device to your head using the velcro straps.
 step 3: Long press the power button for 3s to start the device. It runs automatically for 20 mins.	 step 4: Remove the mask and apply moisturizer to the skin.

Troubleshooting

Indication	Solution
Device does not turn on	*Long press 3 Seconds on ON/OFF switch to turn on *The device needs to be charged
Battery low voltage	The device needs to be charged
USB plug/socket is bad	*Re-insert *Replace the USB plug/socket
Auto power-off	*Long press 3 Seconds on ON/OFF switch to turn on *The device needs to be charged

Operation

on / off	Long press 3 Seconds to turn on/off
Automatic mode	Upon startup, the device enters automatic mode and activates a comprehensive 20-minute light therapy cycle. The cycle consists of three stages: Stage 1 (10 mins) - Infrared (invisible) and red-light therapy for improved circulation and collagen production Stage 2 (5 mins) - Blue light therapy for acne reduction and skin purification Stage 3 (5 mins) - Orange light therapy for increased skin radiance and overall skin health This automated feature allows users to effortlessly enjoy a complete light therapy session in the comfort of their own home.
Manual mode	While in operation, simply select the color and press the "T" key to adjust the working time of the device. With each short press of the key, the working time changes between three options: 15 minutes, 20 minutes, and 30 minutes. This convenient feature allows users to easily customize their light therapy sessions based on their individual needs and preferences.
Adjust Intensity	In working condition , short press "+" "-" key to adjust the brightness (intensity) of the lights.
Charging	When the battery indicator on your device starts blinking, it means that the battery is running low and needs to be charged. Please connect your device to a power source as soon as possible to ensure uninterrupted use.
Use	It is recommended to use skin rejuvenation mask 3-5 times a week

Product specification

Product Name	Beauty LED Mask
Working voltage	3.7V
Charging voltage	DC5V $\pm$ 1A
Rated power	4.5W
LED PCS	368
Product Size	315*219*5mm
Net Weight	$\approx$ 230g
Wavelength	NIR: 850 $\pm$ 5nm, Red light: 630 $\pm$ 10nm, Blue light: 415 $\pm$ 5nm, Orange light: 605 $\pm$ 5nm
Cell Rated Capacity	2000mAh
Material	Food Grade Silicone
Accessories	Face mask, Straps, Controller, Cable, USB cable, Eye protector, User manual

 Please don't use a power charger that is more than 5V, flash charging or fast charging, which will burn out the controller.

People unsuitable to use

1. This product has no fireproof function, please keep away from fire or corrosive liquid;
2. Do not disassemble or repair this product by yourself;
3. Please keep it at a place out of reach of children;
4. This product has waterproof function;
5. Do not press it heavily so as to avoid deformation and failure.

Warranty

The manufacturer guarantees its manufactured products for a period of 24 months from the date of purchase. Accessories included with the original product are warranted for a period of 12 months from the date of purchase.

For the efficient operation of products and complaints, we recommend.

- ① Study and follow the operating instructions carefully.
- ② Keep the manual, warranty card, proof of purchase.
- ③ Use the original accessories supplied with the device.
- ④ In case of problems with the device, contact the Manufacturer.
- ⑤ Do not make repairs or modifications yourself.
- ⑥ Use the product according to its intended use.

During the warranty period, all repairs are performed free of charge.

Any defects or faults will be removed within a maximum period of 14 days, calculated from the time of delivery of the device to the manufacturer's authorized service center.

Limitations

The warranty will not be recognized, in case of mechanical damage due to inattention or fortuitous events.

If the device is used contrary to its intended use and specifications.

Batteries are consumables, the warranty covers 6 months. Batteries cannot be replaced outside of an authorized service center.

Warranty card

Product Name	Beauty LED Mask	Model	Name:	Address:	Tel:	Cause of failure

*Please fill out the warranty card in detail when repairing. Thank you for your cooperation!

Contact details

DIFORO Brand

CUSTOMER SERVICE
+48 453 271 100

Consultation hours : 08:00 ~ 16:00
(Monday to Friday)

PRODUCENT

DIFORO by VOBE sp. z o.o.
Rudna Mała 167A, 36-054 Mrowla, POLAND

www.diforo.pl get@diforo.pl
NIP: 5252799270, KRS 0000800189

The Power of Light and Infrared – Pure Skin Renewal

LED light therapy is a breakthrough in skincare, combining advanced technology with the body's natural regenerative potential. Each wavelength of light penetrates the skin at a specific depth, stimulating its own repair mechanisms – safely, gently, and effectively.

This non-invasive method is widely used in professional cosmetic studios and aesthetic clinics around the world. DIFORO technology now brings this experience to your home – for visible, long-lasting results without irritation or downtime.

Red Light – Revitalization and Youthful Glow

Red light at approximately 630 nm energizes skin cells by stimulating mitochondria to produce more ATP – the “fuel” of cellular regeneration. It encourages collagen and elastin synthesis, improving skin firmness, smoothness, and elasticity. With regular use, fine lines and wrinkles become less visible, circulation improves, and the complexion gains a youthful radiance. Red light is ideal for tired, mature, or dull skin that needs renewal and brightness.

Near-Infrared Light – Deep Repair and Cellular Renewal

Near-infrared light (around 850 nm) reaches the deepest layers of the skin, where it enhances microcirculation, boosts oxygen flow, and accelerates tissue regeneration.

It helps soothe muscle tension, reduce inflammation, and stimulate collagen remodeling for visibly healthier, more resilient skin. This wavelength works from within – awakening the skin's natural ability to repair and rejuvenate itself.

Blue Light – Purity and Balance

Blue light at around 415 nm offers powerful antibacterial properties. It helps neutralize acne-causing bacteria, purify pores, and balance sebum production. It's perfect for oily or combination skin, visibly reducing imperfections and calming redness. The complexion appears cleaner, smoother, and refreshed with every session.

Orange Light – Radiance and Vitality

Orange light (approx. 605 nm) gently energizes dull, fatigued skin, helping to even out tone and restore natural luminosity. It stimulates surface renewal, supports circulation, and reduces signs of stress and tiredness.

Your skin looks brighter, fresher, and full of vitality – glowing naturally, even without makeup.

Synergy of Light – Harmony for Timeless Beauty

By combining multiple wavelengths, Diforo creates a synergistic light experience that works on every layer of the skin. Each session leaves the complexion visibly smoother, healthier, and more radiant.

Regular LED therapy helps maintain youthful firmness, balanced tone, and natural glow – empowering your skin to look its best every day.

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