
Multi functional

Waist massager

Massage belt

Dif-Injury

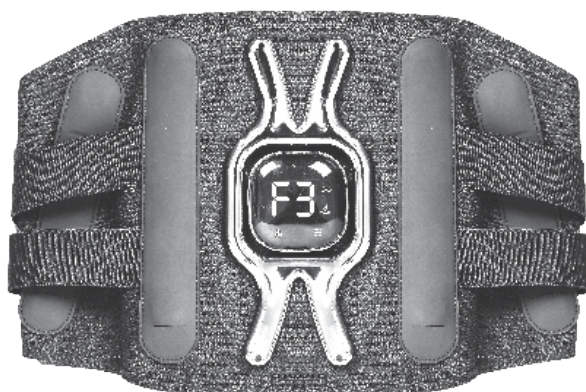operation instructions

1 User must read

Please read and keep this manual carefully before using the product.

Symbol description

To avoid personal injury or property damage, please be sure to follow the following instructions.

- Use the following diagrams to illustrate the desired content to be followed.
-

Indicate prohibition

Indicating the need to follow instructions

Intended for

The following individuals should use it with caution and consult a doctor before use:

- Pregnant women or those who have recently given birth.
- People who are deemed to have lost their sense of warmth.
- Diabetes patients with perceptual impairment.
- People with a body temperature above 38°C (fever period).
- Individuals with skin sensory disorders or skin injuries.
- Those who use medical electronic devices such as pacemakers embedded in their bodies.
- Those who are receiving medical treatment, especially those who feel abnormal physically.
- Patients with lower back injuries and those with lower back sequelae due to traffic accidents, etc.
- Patients with malignant tumors, heart disease, hypertension, and acute illnesses.
- People who have metal, silicon, or other plastic materials implanted in their bodies due to fractures or plastic surgery.
- Osteoporosis patients, those who suffer from severe fractures, sprains, muscle pain, and acute (head) pain.

 The following groups or situations are prohibited from using this product:

- Do not allow individuals with intellectual disabilities or minors to use it alone.
- Minors should use it under the guidance of an adult guardian

Precautions for use

-  It is recommended to use 1-level hot compress and soothing mode at the beginning of using this product.
-  It is recommended to use for 30 minutes each time. Please check your waist condition and adjust the frequency of use appropriately.
-  After using this product, the skin on the waist may appear red for a short period of time, which is a normal phenomenon. If this symptom persistsIf there is no relief, it is recommended to consult a doctor.
-  If you experience any abnormalities or sensations during use, please stop using immediately and consult a doctor.
-  It is recommended not to overuse the massage function.
-  Please use this product in an environment of 5-40 °C.
-  Do not use this product while sleeping.
-  Do not engage in vigorous exercise when using this product.
-  Under no circumstances should you stand or sit on the product.
-  Do not use in a damp environment to avoid electric leakage and injury accidents caused by the product being damp.
-  When using, do not cover it with a blanket or pillow to avoid burns or even fires caused by overheating.


Maintenance precautions

-  When cleaning, be sure to unplug the adapter from the household power outlet.
-  Please handle the packaging bag of this product properly and keep it away from infants and young children.
-  This product contains a lithium battery. Please do not disassemble, modify, immerse the product (controller) in water, or throw it into a fire by yourself, as this may cause damage to the product.
-  Do not store in contact with children and infants. In addition,Do not allow them to come into contact with or operate this product. Otherwise, it may cause burns, electric shock, or personal injury.
-  Do not damage or modify, or forcefully bend or pull the charging cable. Also, do not place heavy objects or squeeze the charging cable on top of it.

2. Product Features


This product uses vibration massage combined with hot compress function to help relieve discomfort in the waist. It has multiple temperatures and massage intensities to choose from, suitable for different people to use for waist massage and relaxation.

3. Component name

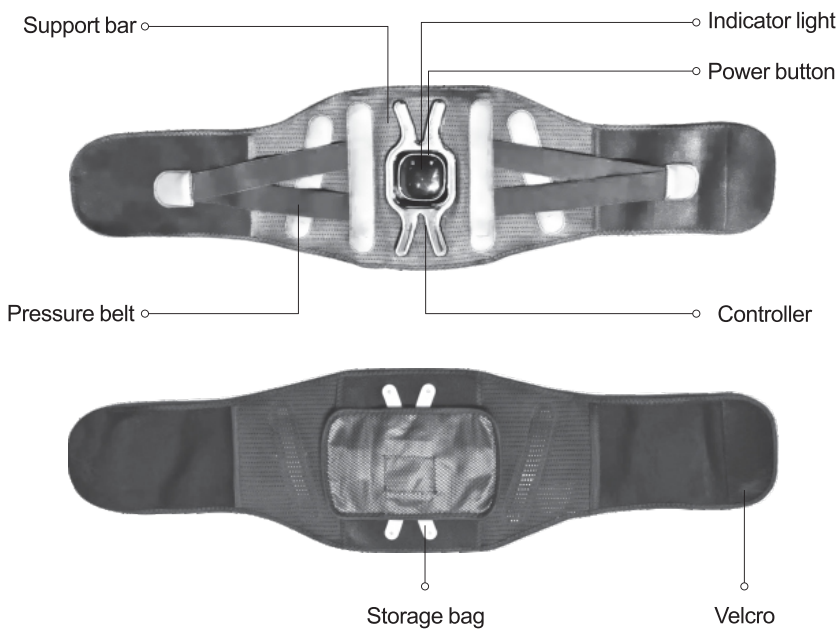
4 Usage Methods

1. Remove accessories such as belts.
2. If you need to use it with a velvet bag, please place the velvet bag in the inner storage bag.

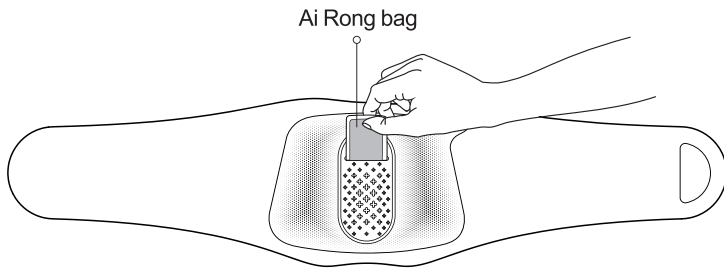2. Adjust the elasticity with Velcro, wear the massager on the waist, and tighten the compression band. Please ensure that the controller is in the lumbar position after wearing it.

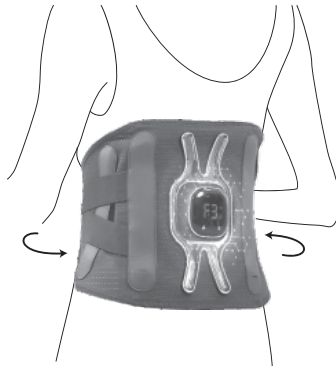- !** Due to the sensitivity of the abdominal skin, it is recommended to use it under clothing when the temperature is 45 °C or above to prevent burns.

Usage operation

Turn on/off

1. Press and hold the power button for 2 seconds to turn on the device, accompanied by a "beep~" long sound prompt. The device defaults to 3 levels of hot compress when turned on.
2. It automatically shuts down within 15 minutes by default, accompanied by a long beep prompt. Switching gears during operation does not affect timing.

Hot compress regulation


Press the heating button to switch the hot compress mode, and cycle through: 1st gear (40°C) -2nd gear (45°C) -3rd gear (50°C). Each operation is accompanied by a "beep" short sound prompt; Display: H1H2H3

Vibration mode adjustment

Press the function keys to switch vibration modes, and cycle through them in sequence: Comfort mode - Rhythm mode - Press mode - Turn off vibration. Each operation is accompanied by two beeps; Display: F1 F2 F3

Charging method

1. Insert the charging cable plug into the charging hole on the side of the controller.
2. Connect the charging cable USB port to a 5V $\overline{=}$ 1A charger and plug it into a power outlet.
3. When charging, turn on the red indicator light, and when fully charged, turn on the green indicator light.
4. It takes about 4 hours to fully charge.
5. Do not turn on the device while charging.


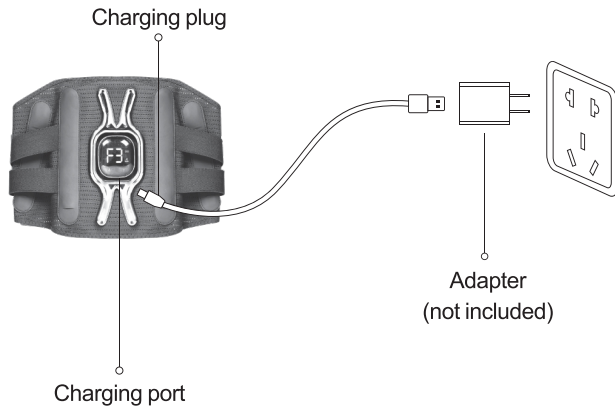
Hot compress gear indicator	State	Describe
H1	Chang Liang	First gear, 40℃
H2	Chang Liang	Second gear, 45℃
H3	Chang Liang	Three gears, 50℃

Vibration mode indication	State	Describe
F1	Chang Liang	Soothing mode
F2	Chang Liang	Rhythmic mode
F3	Chang Liang	Press mode

5 Maintenance and upkeep

Cleaning method

1. Before cleaning, please confirm that the power switch has been turned off and the charging cable plug has been unplugged from the charging port.
2. The control box should be properly protected, and when cleaning the display screen, please use a soft cloth to wipe it;
3. The fabric area supports water flushing. When cleaning and maintenance are needed, it can be wiped or rinsed with clean water, and after rinsing, it needs to be fully dried before use;
4. When washing the fabric area, please remove the massage component first to avoid damage to internal components caused by water ingress;
5. Do not use corrosive cleaning agents such as alcohol, diluents, gasoline, etc. to clean this product, otherwise it may cause the body to fade, and the product structure to corrode and crack.

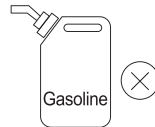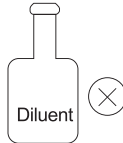When maintaining, be sure to unplug the charging cable, but do not use wet hands to unplug the charging cable.

Disposal method

The lithium-ion battery used in this product is a valuable resource that can be recycled and regenerated. Please remove the lithium-ion battery after use and hand it over to a professional organization that can recycle and regenerate rechargeable batteries for processing.

Battery removal steps:

1. Pry open the top cover of the control box from the bottom;
2. Loosen the screws on the circuit board and remove the circuit board;
3. Cut off the battery connection wire and wrap the wire fracture with insulating tape;

- ❗ Please do not replace the battery by yourself.
- ❗ Before discarding this product, please remove the battery.
- ❗ Please cut off the power before removing the battery.
- ❗ Batteries should be disposed of safely.

6 Fault diagnosis

This chapter lists the most common questions about this product. If you are unable to find the answer to the problem you are facing, please contact the service hotline or authorized service center in your area.

Fault description	Possible reasons	Solution method
Unable to power on	Dead battery	Please recharge
	The product is still connected to the power outlet	Remove the power cord from the body
The charging indicator light does not light up	The USB charging cable is not properly inserted into the power outlet	Please insert the USB charging cable correctly into the power outlet
After charging, it can only be used for a few minutes	Battery life exhausted (approximately 5 years)	Please contact customer service for assistance (please do not replace the battery by yourself)
The body is not heated	When used outside the recommended ambient temperature range, this product may not function properly	Please use in an environment of 5-40 °C
The product stops working during operation	The battery is fully depleted	Low battery, please charge
	30 minute timed end, product automatically shuts down	If you need to continue using, please restart the device

After a long period of disuse, it is necessary to confirm that the function is normal before using it again.		
Regular inspection	Are there any of the following faults	If the above malfunction occurs, in order to prevent accidents, please stop using it immediately and cut off the power supply. Please send it to our company's repair center.
	·Charging cable overheating ·The charging cable has deep scars or deformations, and has a burnt smell ·Abnormal sounds or vibrations during work ·After turning on the switch, it does not work properly ·There are other abnormalities or malfunctions	

be careful:

If the fault cannot be eliminated or the product still cannot operate satisfactorily, contact the after-sales service personnel and do not attempt to disassemble or repair it by yourself.

7 Main Technical Parameters

Product Name	Multi functional waist belt
Product model	Y-10
Working voltage	3.7V ===
Working power	8W
Charging voltage	5V ==
Charging current	1A
Product size	S code 1000x26mm/L code 1100x26mm
Net weight	0.62kg
Battery capacity	2000mAh
Interface type	Type-C
Power supply mode	Built in battery
Date of manufacture	Refer to the product qualification certificate

8 Environmental Protection Instructions

1. This product meets the requirements of the "Management Method for the Restriction of Hazardous Substances in Electrical and Electronic Products";
2. Within the environmentally friendly usage period, consumers can use it with confidence as there will be no harmful substance leakage, precipitation, or other health issues during normal use;
3. The environmental protection period of our company's products is 10 years, and the "environmental protection period" is only valid when the product is used under normal circumstances as described in this manual.

Component Name	Name and content of harmful substances in the product					
	Lead	Lead	Cadmium	Hexavalent chromium	Polybrominated biphenyls	Polybrominated diphenyl ether
Plastic components	○	○	○	○	○	○
Circuit board (electronic) components	×	○	○	○	○	○
Charging line	○	○	○	○	○	○
Hardware components	×	○	○	○	○	○
Battery	○	○	○	○	○	○

This table is prepared in accordance with the provisions of SJ/T 11364.

O: The content of the hazardous substance in all homogeneous materials of the component is below the limit requirements specified in GB/T26572.

X: Indicates that the content of the hazardous substance in at least one homogeneous material of the component exceeds the limit requirements specified in GB/T26572.

(Note: Some of the components marked with "X" contain harmful substances that exceed the standard due to the current limitations of industry technology, which cannot be replaced or reduced temporarily.)

- The following situations will not be eligible for free services

1. Human damage caused by improper use, such as improper power connection, use of inappropriate accessories, improper installation, failure to follow instructions, incorrect use or negligence resulting in damage. Damage caused by transportation and other accidents. Repairs and modifications not approved by our company. Other damages caused by force majeure (such as natural disasters, abnormal voltage, etc.).

2. Damage caused by non household use (such as industrial or commercial use).

3. Product aging, wear and tear caused by normal use, but not affecting the normal use of the product.

- This certificate is only valid within the territory of China.

·The warranty certificate, as a guarantee of free maintenance within the prescribed period and conditions, does not limit the legal rights of consumers.

|

People unsuitable to use

- 1. This product has no fireproof function, please keep away from fire or corrosive liquid;
- 2. Do not disassemble or repair this product by yourself;
- 3. Please keep it at a place out of reach of children;
- 4. This product has waterproof function;
- 5. Do not press it heavily so as to avoid deformation and failure.

Warranty

The manufacturer guarantees its manufactured products for a period of 24 months from the date of purchase. Accessories included with the original product are warranted for a period of 12 months from the date of purchase.

For the efficient operation of products and complaints, we recommend.

- ① Study and follow the operating instructions carefully.
- ② Keep the manual, warranty card, proof of purchase.
- ③ Use the original accessories supplied with the device.
- ④ In case of problems with the device, contact the Manufacturer.
- ⑤ Do not make repairs or modifications yourself.
- ⑥ Use the product according to its intended use.

During the warranty period, all repairs are performed free of charge.

Any defects or faults will be removed within a maximum period of 14 days, calculated from the time of delivery of the device to the manufacturer's authorized service center.

Limitations

The warranty will not be recognized, in case of mechanical damage due to inattention or fortuitous events.

If the device is used contrary to its intended use and specifications.

Batteries are consumables, the warranty covers 6 months. Batteries cannot be replaced outside of an authorized service center.

Contact details

DIFORO Brand

CUSTOMER SERVICE
+48 453 271 100

Consultation hours : 08:00 ~ 16:00
(Monday to Friday)

PRODUCENT

DIFORO by VOBÉ sp. z o.o.
Rudna Mała 167A, 36-054 Mrowla, POLAND

www.diforo.pl get@diforo.pl
NIP: 5252799270, KRS 0000800189

Warranty card

	Model				
Beauty LED Mask					
Product Name		Name:			
Distributor information		Address:			
		Tel:			
		Buyer information			Cause of failure

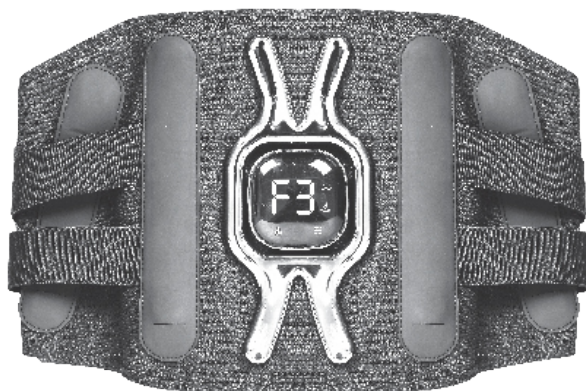
*Please fill out the warranty card in detail when repairing. Thank you for your cooperation!

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