
GIO-460

innoGIO

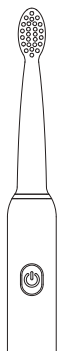GIOsonic
giraffe/ crocodile/ panda/
Coco the capybara

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@innogio_com

Note:


Please read all instructions carefully before using the product.
Please report any errors to us.

GIOsonic
giraffe/ crocodile/ panda/
Coco the capybara

Sonic toothbrush for children

USER MANUAL

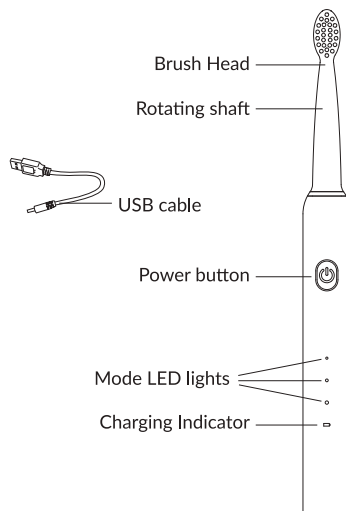Specifications & Information

Item Name: InnoGIO GIOsonic Giraffe/Crocodile/ Panda/ Coco the capybara – Sonic toothbrush for children
Model: GIO-460GIRAFFE/ GIO-460CROCODILE/ GIO-460PANDA/ GIO460-CAPYBARA
Waterproof level: IPX7
Rated Voltage/Power: DC 3.7V–0.5W
Charging Time: about 3 hours
What's in the Box: Toothbrush Handle and Head, USB Cable, User Manual

Dental Hygiene Advice

*Brush your teeth at least twice a day.
Change the brush head every 2-3 months.*

Low level (Primary teeth period)

Cultivate the baby's toothbrushing habits; clean its mouth gently and use soft bristles.

Medium level (Teeth eruption period)

To protect healthy growth of young teeth, take care of sensitive gums and ensure teeth are brushed every day.


High level (Permanent teeth period)

Maintain good toothbrushing habits; clean teeth thoroughly to reduce the risk of tooth decay.


Charging Instructions

1. Check the charger and the socket – make sure they are safe, standard-compliant and appropriate to the voltage used.
2. Make sure the USB cable is dry and the charger plug free from water stains, damage, etc. Plug in the charger in safe conditions.
3. Before using the electric toothbrush for the first time, turn it on. Fully charge the battery before the first use in case of low or discharged battery.


Steps

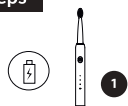When the battery is low and the LED charging indicator is flashing, charge the toothbrush before using it. For information on charging, see the charging instructions section.

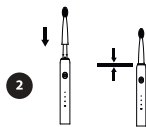
Attach the brush head, as shown in the picture. The gap between the brush head and the handle is normal.

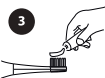
Wet the brush head and apply an appropriate amount of toothpaste on the bristles.


Put the brush head in your mouth; hold the bristles at an 45 degree angle towards the gum line: press the button and select the desired brushing mode.

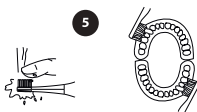Move the brush head up and down; do not press hard or scrub.

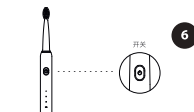To prevent splashing of the toothpaste foam, turn off the toothbrush before taking it out of your mouth.

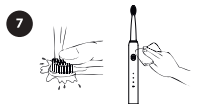Rinse the brush head once you have finished brushing.

Keep the brush head in a dry and ventilated place.

Battery Indicator Description

The light status when not charging	Description
Indicator	Standby Mode
Red light on	Low battery, charging required
Light off	Battery discharged, the device cannot work

Light status when charging	Description
Red light on	Charging Mode
Red light turned green	Charging finished

DISPOSAL OF UNNECESSARY EQUIPMENT BY HOME USERS IN THE EUROPEAN UNION

The presence of this symbol on the product or its packaging indicates that this product must not be disposed of in the same way as household waste. As a result, you are responsible for the proper disposal of the equipment, and you are obliged to deliver it to an authorized recycling point for waste electronic and electrical equipment. Sorting, disposing and recycling your waste equipment will help to conserve natural resources and ensure that it is recycled in a manner that respects the environment and human health. For more information on collection points for used equipment, please contact your local household waste disposal or local authority.

innoGIO

www.innogio.com

Importer: **INNOGIO GLOBAL** sp. z o.o.
Ostródzka 74H, 03-289 Warsaw, Poland

@innoGIO.official

@innogio_com