
medivon.

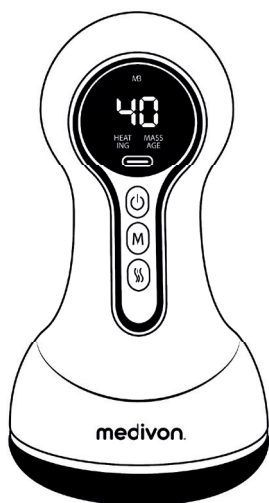User manual

Medivon Bian

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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Device operation

Medivon Bian is a modern massager made of natural stone, combining shockwave therapy, vibration, and heating functions. Since the device is made of natural stone, slight discolorations, lines, or patterns may appear on its surface. These are natural characteristics of the material and do not affect the performance of the massager. The device is fragile – avoid dropping or striking it.

Before first use, it is recommended to wipe the device with a dry cloth, as it may be covered with a protective layer. Before the massage, it is advisable to apply massage oil or lotion to the skin. Start with the lowest intensity level and gradually increase the power over the following days. The optimal massage time is 10–20 minutes per day.

The device is turned on by pressing the power button once and turned off by pressing it again. The massager automatically switches off after approximately 30 minutes of operation. It features 5 massage intensity levels and 3 heating levels, reaching a maximum temperature of about 50°C. Before starting the massage, preheating the device for 5–10 minutes is recommended. If you have sensitive skin, avoid excessive heat.

A full charge takes approximately 4–5 hours, and operating time on one charge is 2–3 hours. During charging, the symbol “--” flashes on the display. When a steady “88” appears, the device is fully charged.

The massager is intended for use on areas such as the neck, shoulders, abdomen, thighs, and legs. It supports muscle relaxation, improved circulation, and overall body recovery.

The product is not waterproof and should not be used in a humid environment or with wet hands. Do not use chemical cleaning agents or disassemble the device. After each use, gently wipe it with a soft, slightly damp cloth.

This product does not replace medical devices or pharmacological treatment. In case of health issues or skin conditions, consult a doctor.

Keep the device out of reach of children. The device is not a toy. Individuals with limited physical or cognitive abilities may use the device only under the supervision of another person. Do not modify the device or attempt to repair it yourself. This may result in damage to the device and void the warranty.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 3 August 2005.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the Medivon Bian product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU as regards the list of restricted substances.

The full text of the Declaration of Conformity is available at the following website:

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Importer: Armare S.A.
ul Grzybowska 87,
00-844 Warszawa, Poland
Made in China
Marka: Medivon
Model: Medivon Bian
E-mail: serwis@medivon.pl
www.medivon.pl
SN:01265904119283550

Producent: Wenzhou Bolang Beauty & Hairdressin
Equipment Co., Ltd, Xihu Village, Nanbin Street,
Rui'an City, Wenzhou City, Zhejiang Province, China
E-mail: selina@micro-beauty.com.cn

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>