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Oriona

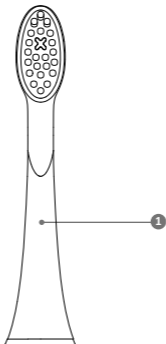
Table of contents


Spis treści

5	Instrukcja obsługi	PL
6	User manual	EN
8	Bedienungsanleitung	DE
10	Návod k použití	CZ
12	Használati utasítás	HU
13	Používateľská príručka	SK
15	Användarhandbok	SE
17	Käyttäjän käsikirja	FI
18	Brukerhåndbok	NO
20	Brugermanual	DK
22	Gebruikershandleiding	NL
23	Manual del usuario	ES
25	Manuale utente	IT
27	Manuel de l'utilisateur	FR
29	Manual de utilizare	RO
30	Korisnički priručnik	HR/BA
32	Uputstvo za upotrebu	RS/ME
34	Uporabniški pravilnik	SI
35	Εγχειρίδιο χρήστη	GR
37	Lietotāja rokasgrāmata	LV
39	Vartotojo vadovas	LT
40	Kasutusjuhend	ET
42	Посібник користувача	UA
44	Ръководство за потребителя	BG
46	Упатство за употреба	MK
48	Manual do utilizador	PT


1


2

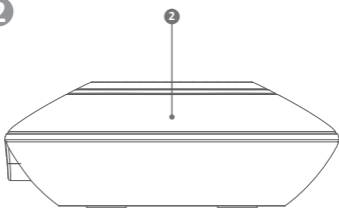
3

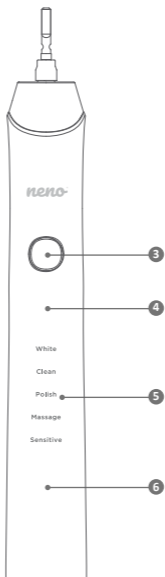4

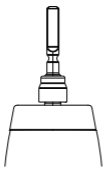5

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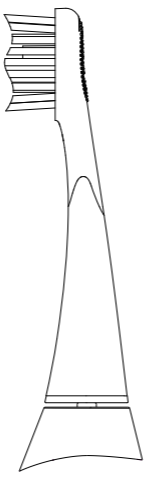EN

USER MANUAL

Dear Customer, thank you for purchasing Neno Orion. The product you have purchased is an electric toothbrush for adults. This device is used to clean teeth and maintain oral hygiene. Please read the following instructions before using the device.

01. PRECAUTIONS

1. Do not disassemble the product. This may damage the product, cause electric shock and invalidate the warranty.
2. Be careful not to drop the device on the floor. If it has been damaged, contact the manufacturer.
3. Do not use accessories other than those recommended by the manufacturer.
4. Do not use flammable or explosive substances to clean the product.
5. This product is intended for brushing teeth only. Do not use the device for any other purpose.
6. Discontinue use if the device is visibly damaged.

7. Do not place the product in water. This may cause malfunction and pose a risk of fire, electric shock or explosion.
8. Keep the device away from water containers such as bathtubs or sinks.
9. If the toothpaste used contains peroxides or bicarbonates (often found in whitening toothpastes), wash the device thoroughly after use to avoid damaging its plastic casing.
10. The device is not intended for use by persons with limited physical, cognitive or mental abilities or who lack experience or knowledge. If such a person wishes to use the device, they must be supervised or instructed on how to use it by a responsible person.
11. The device is not a toy.
12. Do not expose the toothbrush to direct sunlight.
13. The device is for personal hygiene and should not be used by more than one person.
14. Do not use the charging stand outdoors or in places with high temperatures.
15. Do not expose the product to direct sunlight or place it near heat sources.
16. Do not immerse the base in water or any other liquids.
17. Before using or cleaning the product, disconnect it from the power source.

02. DESCRIPTION OF THE DEVICE

Figure 1

1. Brush head

Figure 2

2. Charging stand

Figure 3

3. Mode selection button
4. Toothbrush handle
5. Mode indicator
6. Charging indicator

03. SET CONTENTS

1. Toothbrush grip
2. Two brush heads
3. Charging stand
4. User manual
5. Travel case

04. HOW TO USE

1. Place the brush head (Figure 1-1) on the handle (Figure 3-4) as shown in Figure 4.
2. Wet the brush head (Figure 1-1) with water and apply toothpaste.
3. Gently place the bristles of the brush against your teeth, along the gum line.
4. Select the appropriate mode by pressing the mode selection button (Figure 3-3). The mode indicator (Figure 3-5) will show which mode has been selected. Every 30 seconds, the toothbrush will pause for 0.5 seconds to remind the user to change the brushing area. The total brushing time is 2 minutes — after this time, the toothbrush will automatically switch off.
NOTE: You can select one of the following modes:
 - a) **White** – The whitening mode removes surface stains, brightening your teeth.
 - b) **Clean** – A very gentle mode for sensitive teeth and gums. Recommended for daily use.
 - c) **Polish** – A quick 2-minute mode that whitens and polishes your teeth.
 - d) **Massage** – The massage mode gently stimulates and massages the gums, contributing to their better condition.
 - e) **Sensitive** – Some people may experience discomfort when using it for the first time. The vibrations of the bristles vary from weaker to stronger to make it easier to get used to using a sonic toothbrush.
5. After brushing your teeth, rinse the toothbrush under running water. Shake off any remaining water from the bristles and use a soft cloth to wipe the handle of the device.

05. CHARGING THE TOOTHBRUSH

1. Place the toothbrush (Figure 3-4) on the charging base (Figure 2-2), which is connected to a power source. The toothbrush will begin charging via inductive charging.
2. Charging is in progress when the LED is lit (Figure 3-6). When fully charged, the LED turns off. A full charge takes no longer than 24 hours.
3. The mode selection button (Figure 3-3) does not work during charging.
4. If the brush handle is in operation and is placed on the charging base, the brush will stop operating and enter charging mode.

06. CLEANING THE TOOTHBRUSH

1. **Handle:**
 - Remove the brush head (Figure 1-1) from the handle (Figure 3-4) and wipe it with a soft, damp cloth, then

leave it to dry.

- Do not immerse the handle in water for long periods of time.
- 2. Charging base:**
 - Wipe the charging base (Figure 2-2) with a soft, damp cloth and leave it to dry.
 - Do not immerse it in water.
 - 3. Brush head:**
 - After each use, rinse the bristles and head (Figure 1-1) with water, then leave to dry.
 - Do not bend the bristles forcefully to prevent them from losing their shape.
 - It is recommended to replace the brush head every three months. After a period of use, the bristles may become deformed, rough and prone to accumulating bacteria, which affects cleaning performance.

07. SPECIFICATION

Power supply: 5V – 0,3A

Battery: Lithium-ion 700 mAh

Water resistance: IPX7

Brush bristles: Nylon

Vibrations per minute: up to 31,000

Built-in timer: yes

Dimensions: 2.7 x 2.7 x 25.4 cm

Weight: 117g

Full charging time: 24h

08. WARRANTY CARD

The product is covered by a 24-month warranty. The warranty terms and conditions can be found at: <https://nenopl/gwarancja>

Details, contact information and service address can be found at:

<https://nenopl/kontakt>

Specifications and kit contents are subject to change without notice.

We apologise for any inconvenience.

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EN

The crossed-out wheeled bin symbol indicates that unwanted electrical or electronic equipment, including accessories (such as power supplies or cables) and components (e.g., batteries, if included), must not be disposed of with household waste. Instead, these items should be taken to an appropriate collection point where they will be accepted free of charge for proper disposal or recycling. This requirement is in accordance with the revised WEEE Directive (2012/19/EU) and Regulation (EU) 2023/1542, which amends Directive 2008/98/EC and Regulation (EU) 2019/1020, and repeals Directive 2006/66/EC. Proper disposal of such devices helps prevent environmental damage. Information on local collection points is available from the relevant municipal authorities. Improper disposal of waste may be subject to penalties under applicable local laws.

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