TENS/EMS PHYSIOTHERAPY MASSAGER

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1. GENERAL INFORMATION

1.1 What is this device?

This portable electrostimulator combines TENS, EMS, and massage functions. Before first use, read the entire instruction manual carefully and retain it for future reference.

1.2 A few words about pain

Pain acts as an alarm in our bodies – it signals a problem, helping us avoid more serious bodily harm. Without it, minor abnormalities could go unnoticed and lead to injury. However, chronic pain loses its usefulness and becomes a source of suffering. The pain signal originates at the site of injury, travels through small nerves to the spinal cord, and from there travels to the brain. There, it is processed, interpreted, and experienced as discomfort.

1.3 What is TENS?

TENS, or transcutaneous electrical neural stimulation, is a proven method of pain relief. It's used daily by physiotherapists, medical personnel, and professional athletes worldwide.

High-frequency impulses block the transmission of pain signals to the brain, activating the nervous system's natural defense mechanisms. Electrodes placed on the skin in the area of pain stimulate the nerves, which suppresses the sensation of pain. Low-frequency TENS currents also stimulate the release of endorphins—the body's own analgesics.

1.4 What is EMS?

EMS (electrical muscle stimulation) is a globally recognized technique for rehabilitating muscle injuries. It sends electrical impulses directly to muscles, forcing them to contract without any active effort.

The technology is based on the square wave, invented in 1831 by John Faraday. The low frequency and waveform allow for precise action on motor neurons, engaging entire muscle groups.

1.5 What is the RELAX function?

Massage is not a medical treatment, but a relaxing vibration therapy. This program gently loosens tense muscles, providing relief and improved comfort.

2. SAFETY INFORMATION

2.1 Application of the device

TENS mode

It is used for short-term pain relief in the muscles of the neck, shoulders, back, joints, hips, hands, abdomen, feet, and upper and lower limbs. It also helps with pain resulting from physical activity or daily household chores. It also supports the treatment of chronic, difficult-to-control pain and alleviates symptoms associated with arthritis.

EMS mode

- Maintaining or expanding range of motion
- Relaxing tight muscles
- Preventing or delaying muscle atrophy

- Retraining muscle function
- Improving local blood circulation
- Rapid post-operative calf muscle stimulation to prevent venous thrombosis

RELAX mode

Contains 12 programs to relax and support muscle building

Important! The MT6543 is NOT a medical device and cannot replace treatment prescribed by a qualified physician. The device is intended for home use only, by adult users.

Important! If you suffer from or suspect any medical conditions, consult your doctor before using MT6543.

3. WARNINGS AND SAFETY

3.1 Key Warnings and Safety Precautions

Before using this device, please read the following guidelines carefully. They are intended to protect your health, minimize the risk of injury, and prevent equipment failure.

3.2 Prohibited

- Do not use on pacemakers, defibrillators, or other metal/electronic implants – risk of shock, burns, interference, or death.
- Avoid areas with tumors or other skin lesions.
- Do not use on swelling, inflammation, or eruptions (e.g., varicose veins, thrombosis).
- Electrodes outside the carotid sinus and over the head.
- Do not use on hyperactive tissue, inguinal hernias, or postoperative scars (min. 10 months old).
- Do not use on severe circulatory disorders in the legs.

3.3 Pay close attention

- Consult a doctor before use if the pain is being treated with medications or therapies.
- Discontinue use if the pain worsens, does not subside after 5 days, or is severe – consult a doctor.
- Do not stimulate the neck (risk of cramps, shortness of breath, arrhythmia).
- Not on the chest (risk of fatal arrhythmia).
- Avoid tumors, medical equipment (ECG), bathtubs/showers, while sleeping, or driving.
- Only use on healthy skin; does not replace medication.
- Not near the heart, eyes, face, genitals, or insensitive areas; keep electrodes separate.
- Store out of reach of children; consult a doctor if you experience any concerns.
- Do not increase intensity if discomfort occurs.


3.4 Be careful

- TENS does not treat central pain (e.g., migraines) or replace medication – it is a symptomatic therapy.
- Safety during pregnancy has not been confirmed; avoid overuse.
- Risk of skin irritation from the gel/stimulation; use with caution in cases of epilepsy, heart disease, bleeding, or post-operative conditions.
- For single-user use only; not for people with dementia/low IQ.
- Check electrodes; use only recommended ones; not with other electrical impulses.
- Do not use sharp instruments on buttons; not during menstruation/pregnancy without consulting a doctor.


3.5 Adverse effects

- Possible irritation/burns under the electrodes.
- Rare: dizziness/fainting in novices – start sitting.
- Reduce intensity if discomfort occurs; consult a physician if it persists.

4. ABOUT THE PRODUCT

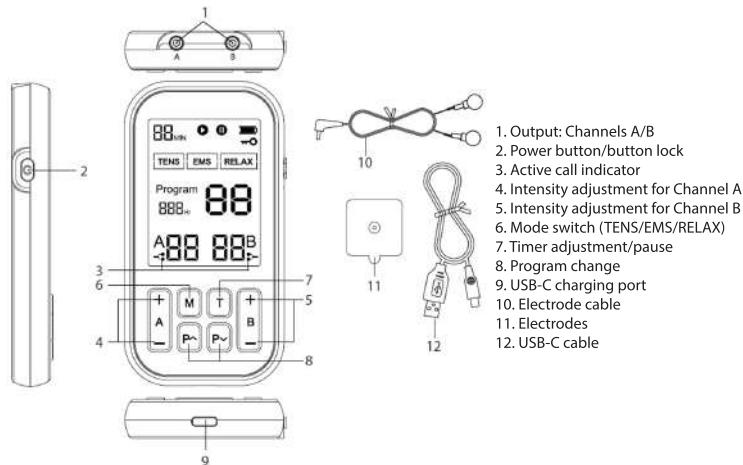
5. LCD DISPLAY

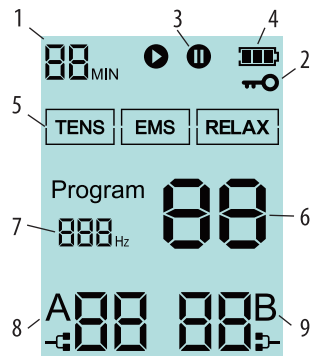- Zone 1: Treatment Time: Displays the remaining treatment time.
- Zone 2: Button Lock Indication
- Zone 3: Pause Indication during treatment
- Zone 4: Battery Level: When turned on, the battery indicator on the outer frame lights steadily, and the internal three-segment battery level display shows the battery level. As the battery level drops, the corresponding number of segments decreases sequentially from right to left. When the battery level reaches a critical level, a low battery alarm appears, the outer frame flashes for 10 seconds, and the product automatically turns off.
- Zone 5: TENS/EMS/RELAX Modes. The default mode upon startup is TENS mode. The mode icon flashes to indicate the selected large mode of the current channel.
- Zone 6: Current Mode Program Number Display Area: The default number at startup is 1. TENS mode has 24 programs. EMS mode has 24 programs. There are 12 programs available in RELAX mode.
- Zone 7: Frequency display area: The value reflects the currently generated pulse frequency from the channel.
- Zone 8/9: Displays activity for channels A and B.

6. USE


1. Connect the electrode leads to the A and B channel sockets (if you only want to stimulate one point, connecting to channel A is sufficient)
2. Connect the electrodes to the snaps at the ends of the cables. Peel off the protective film, being careful not to touch the adhesive part with your fingers.

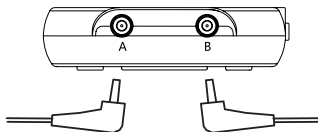Important! The device must be completely OFF.


Routinely replace electrodes when they lose their ability to adhere properly. If you have any concerns about their adhesion, order new replacement electrodes immediately.

Replace them according to your physician's or device manufacturer's instructions to ensure high-quality stimulation. When applying fresh electrodes, carefully follow the application instructions on the packaging – this will ensure effective therapy and avoid skin irritation.

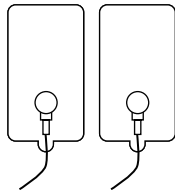3. Place the electrodes on the areas to be stimulated, following the instructions in this manual. It's also helpful to refer to the acupuncture point guide.


The areas to which the electrodes will be attached must be clean, dry, and as hairless as possible.

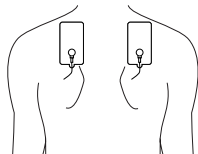
Important!

- **Never place electrodes on wounds or areas with irritated skin.**
- **Remove electrodes from the body with a gentle, even motion, especially on very sensitive skin, to prevent injury.**
- **To prevent possible skin irritation, apply a soothing cream to the electrode areas.**
- **Never turn on the device if the adhesive electrodes are not properly secured in the clasps.**
- **To remove or reposition the electrodes, first turn off the entire device or the selected channel – this will prevent unwanted irritation.**
- **Adjust the electrode size to the size of the area being stimulated.**
- **Never remove the adhesive electrodes while the device is in operation.**


4. Turn on the device by pressing the power button

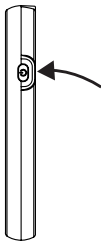
5. Press „M“ to select the desired operating mode

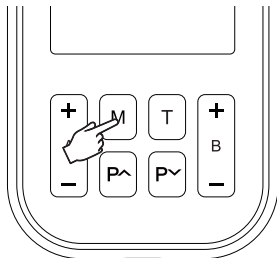6. Use the „P” buttons to select the therapy program

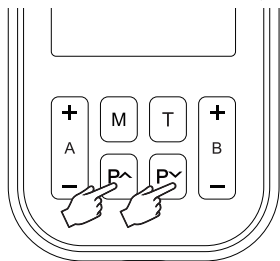7. Use the +/- buttons to adjust the pulse intensity for channels A and B.

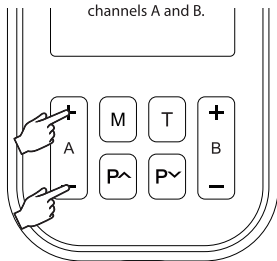Important! The intensity of the pulses must be set so that they do not cause pain, but only a tingling sensation in the treated area.

8. Short presses of the „T” button set the treatment time. Long presses pause or resume the treatment.

7. MAINTENANCE AND CARE

1. After using the device, turn it off and remove the electrodes.
2. Applying a few drops of clean water before and after use will extend the life of the electrodes. Do not use any cleaning agents or detergents.
3. Stick the electrodes back onto the foil from which they were removed to prevent them from drying out. We recommend storing them in zip-lock bags.
4. Replace the electrodes after 15-20 uses (this depends largely on use, body type, and care).
5. For hygiene reasons, each user should use their own set of electrodes.
6. Store the device in a cool, dry place until next use.
7. Clean the device with a clean, slightly damp cloth without detergents.
8. Turn off the device when cleaning.
9. Charge the device at least every 3 months if not used regularly.

8. TECHNICAL DATA

- 3-in-1 physiotherapy device using electrical impulses
- Relaxes, helps relieve pain, and promotes blood circulation
- Easy and safe to use
- Pocket-sized, large, easy-to-read LCD screen
- Intuitive button control
- 60 programs:
 - TENS – nerve stimulation and pain relief – 24 levels
 - EMS – muscle stimulation – 24 levels
 - Massage mode – 12 levels
- Treatment time adjustment: 10-60 minutes
- 2 channels A/B, 2 electrodes per channel – simultaneous therapy in two locations
- 3 sets of electrodes, with self-adhesive hypoallergenic gel
- Medical-grade gel, latex-free
- Power supply: 180mAh Li-ion battery, up to 60 minutes of operation
- Automatic shut-off when inactive
- Charging: USB-C port
- Dimensions: 110x60x15mm
- Weight: 69g

8. PROBLEM SOLVING

Problem	Possible cause	Solution
One electrode seems stronger than the other.	This is normal. Different areas of the body react differently.	There's nothing you need to do. Just make sure the electrodes are moist and sticking well.
The intensity seems very weak, even on the highest settings.	The electrodes are not securely attached to the body. The clear films are still on the electrodes. The lead is not properly connected to the device.	Firmly attach both electrodes to the skin. Remove the protective film from the adhesive surface of the electrodes. Do not allow the electrodes to overlap. Securely plug the cable into the socket.
The intensity seems to decrease over time at the same level	Battery capacity is low.	Increase the intensity level when the sensation becomes less effective. Recharge the battery.
The skin turns red or is irritated	The adhesive surface of the electrodes is dirty or dry. The therapy duration is too long or the intensity is too high. The electrode surface is worn.	Gently wash the adhesive surface of the electrodes with your fingertips for approximately 3 seconds under slowly running water. Reduce the duration of use or the intensity. Replace the electrodes.
No power; nothing displayed on the LCD	The battery is discharged.	Charge the battery.
Power turns off during use.	The battery is dead. The cable is damaged.	Charge the battery. Replace the cable.
It is difficult to attach the electrode to the skin	The protective film on the adhesive surface of the electrodes was not removed. The electrode was applied to wet skin. The adhesive surface of the electrode is damaged.	Remove the protective film from the adhesive surface of the electrodes. Dry the electrode. Replace the electrode.
The adhesive surface of the electrode does not stick	The electrodes were damaged. The electrodes were stored in high temperatures, humidity, or direct sunlight.	Contact us for a replacement. Replace the electrode.


9. PLACEMENT OF ELECTRODES

Below we have included recommendations for electrode placement for each mode.

Important! The electrode size should be matched to the area being treated.


Important! The recommended time for each treatment (per area) is 15 minutes. Do not exceed 30 minutes.

TENS MODE – PAIN REDUCTION

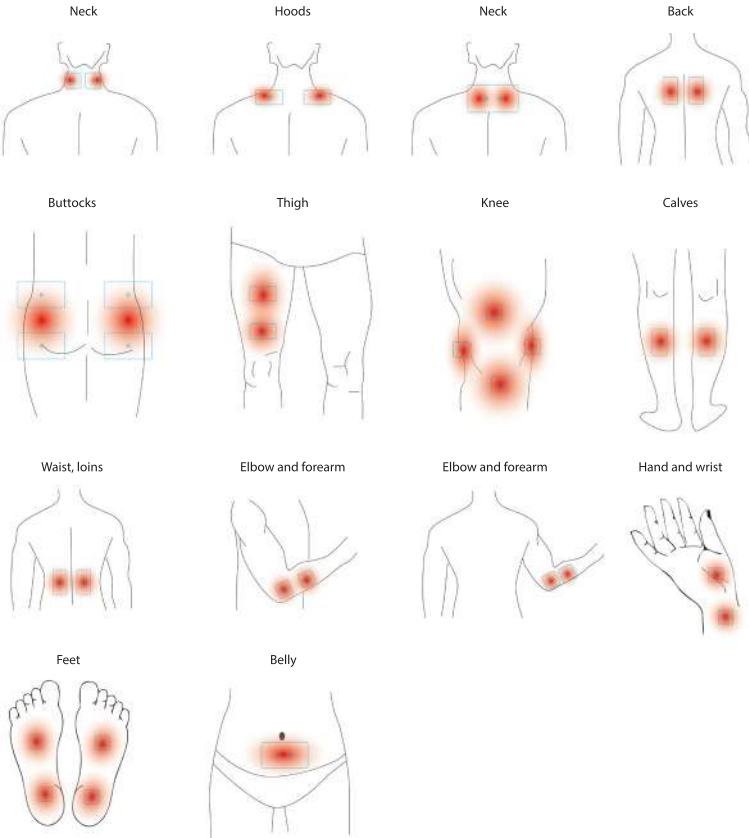TENS MODE – PAIN REDUCTION

Back – trapezoidal

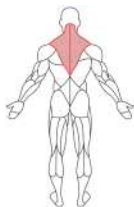
Front – shoulder straps

Center – deltoid


Back – shoulder

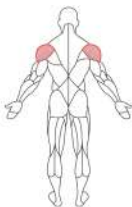
Front – biceps


Back – triceps

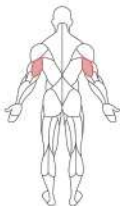Back – latissimus dorsi

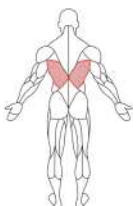Back – trapezius lumbar

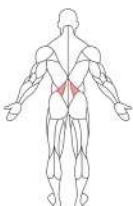
Front – rectus abdominis


Front – internal and external abdominal obliques


Back: gluteal muscles

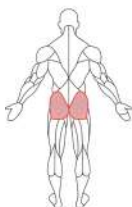Back: gluteal muscles

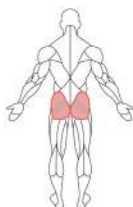Back – large and medium gluteal muscles

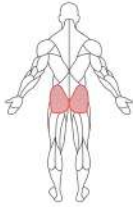Front – rectus femoris, quadriceps femoris, medial femoris, lateral femoris

Front – rectus femoris, quadriceps femoris, medial femoris, lateral femoris


Back – ischio-tibial (biceps) of the thigh

Front – tibialis anterior

Back – Gastrocnemius

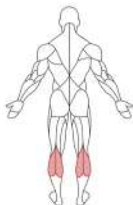RELAX MODE – MUSCLE RELAXATION

In this mode, we recommend using the TENS and EMS point map. The pulses are designed to relax specific muscle groups by simulating tapping and pulsing.

In all modes, we recommend using the acupuncture point map.

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EN

The device is marked with a crossed-out wheeled bin symbol, in accordance with European Directive 2012/19/EC on waste electrical and electronic equipment. Products marked with this symbol should not be disposed of or thrown away from other household waste at the end of their useful life. The user is obliged to dispose of waste electrical and electronic equipment by delivering it to a designated point where such hazardous waste is recycled. Collecting this type of waste in separate locations and recycling it properly contributes to the protection of natural resources.

If the product contains batteries covered by Regulation (EU) 2032/1542 of the European Parliament and of the Council, it cannot be disposed of with normal household waste. Please refer to local regulations regarding the separate collection of batteries, as proper disposal will prevent potential negative effects on the environment and human health.

Proper recycling of waste electrical and electronic equipment has a beneficial effect on human health and the environment. For information on where and how to dispose of waste electrical and electronic equipment in an environmentally safe manner, users should contact their local authority, waste collection point, or the point of sale where the equipment was purchased.

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EN - General safety rules

- Before first use and operation of the Media-tech device, read the user manual carefully and follow the precautions below
- The product is not intended for children aged 0-3 years
- The product may contain packaging elements (foils, blisters) that, if used incorrectly, may cause serious bodily harm or even death
- Do not immerse the product in water unless it is waterproof (marked at least IPX7)
- Do not use in a humid environment, unless it is waterproof at IPX4
- Unless otherwise marked, the product is intended for home use at room temperature
- Use only the original accessories supplied in the set. Non-original ones may cause damage or other serious events
- Before starting the device, make sure that all elements are properly installed and secured.
- For products powered by a battery or accumulator: use only batteries that meet the technical requirements of the device
- In the event of damage to the power cable or other electrical part, do not attempt to repair the device yourself.
- Connect the device only to a suitable power source with parameters consistent with the product specifications.
- Unauthorized disassembly of an electrical, electronic or battery-based product, modifications or attempted repair may result in electric shock and/or ignition and fire.
- If an electrical or battery-based product begins to emit smoke, a burning smell, or becomes abnormally hot during operation, immediately stop using it, disconnect it from the power source and put it in a safe place away from flammable materials.
- A device equipped with an O3 ozone generator emits a minimum concentration of this gas – lethal to bacteria and fungi - not dangerous to people and animals. However, we do not recommend direct inhalation directly above the device.
- A device equipped with a UV diode emits a minimum value of these rays - lethal to bacteria and fungi - not dangerous to people and animals. However, we do not recommend observing this light from a very close distance
- In the case of speakers and headphones, too high a volume setting can damage hearing
- If the product contains colored and pulsating lighting, it should be used with great caution in the vicinity of people with epilepsy
- None of the Media-tech products are medical products and may be used as part of medical therapy
- Hygiene products (irrigators, toothbrushes) should be used strictly for their intended purpose in the areas of the body designated for this purpose. Other use may result in serious bodily harm or disability
- The use of mechanical massage products should be subject to medical consultation in the case of people with cardiological, respiratory and bone diseases
- The use of EMS microcurrent massage products should be subject to medical consultation in the case of people with cardiological and skin diseases
- EMS massagers used contrary to the instructions or for too long (above the permissible duty cycles) may cause local, superficial damage to the epidermis
- Pregnant women should consult a doctor about the possibility of using all massage devices

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