MESH NEBULIZER MIST MANUAL

Shanghai International Holding
Corp.GmbH (Europe)
Eiffestrasse 80, 20537 Hamburg, Germany

Dongguan SIMZO Electronic Technology Co.,Ltd.
NO.6, Zhangzhou Road, Daojiao Town, Dongguan City,
Guangdong Province, 523187, P.R.China.

ECN005

MESH NEBULIZER

The aerosol is generated by pushing liquid medicine through fine holes in the newly developed metal alloy mesh. This is realised by a titanium vibrator, oscillating at high frequency.

GENERAL SAFETY INSTRUCTIONS

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible, the box with the internal packing. If you give this device to other people, please also pass on the operating instructions.

- The appliance is designed exclusively for private use and for the envisaged purpose. This appliance is not fit for commercial use.
- Do not use it outdoors. Keep it away from sources of heat, direct sunlight, humidity (never dip it into any liquid) and sharp edges. Do not use the appliance with wet hands. If the appliance is humid or wet, unplug it immediately.
- When cleaning or putting it away, switch off the appliance and always pull out the plug from the socket (pull the plug itself, not the lead) if the appliance is not being used and remove the attached accessories.
- Do not operate the device without supervision. If you leave the room you should always turn the device off. Remove the plug from the socket.
- The device and the mains lead have to be checked regularly for signs of damage. If damage is found the device must not be used.
- Do not try to repair the appliance on your own. Always contact an authorized technician. To avoid the exposure to danger, always have faulty cable be replaced only by the manufacturer, by our customer service or by qualified person and with a cable of the same type.
- Use only original spare parts.
- Pay careful attention to the following "Special Safety Instructions".

CHILDREN AND FRAIL INDIVIDUALS

- In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

WARNING!

Do not allow small children to play with the foil as there is a danger of suffocation!

- This device is not intended to be used by individuals (including children) who have restricted physical, sensory or mental abilities and/or insufficient knowledge and/or experience, unless they are supervised by an individual who is responsible for their safety or have received instructions on how to use the device.
- Children should be supervised at all times in order to ensure that they do not play with the device.

SPECIAL SAFETY INSTRUCTIONS

- Only use the type and amount of medicine as prescribed and instructed by your doctor.
- In some positions (e.g. upright, nebulization may stop after a short while. In this case, briefly tilt the unit again in order to re-immers the vibrator in the medicine.
- Do NOT poke the mesh with a cotton swap or a pin; the mesh may be permanently damaged.
- If the device is not cleaned and disinfected correctly and as frequently as indicated, microorganisms may remain in the unit, causing risk of infection.
- If the device keep working for continuous 20 minutes, it will switch off automatically after

3 times blinks in red. Press the "on/off" button to switch on again.

• Do NOT immerse any other parts of the nebulization unit in water except the inside of medicine container.

SPECIFICATIONS

Power Supply: 1.5V *2AA Alkaline battery (not included) / USB cable

Power consumption: Approx. 1.5 W

Ultrasonic frequency: 103 KHZ~123 KHZ

Standby current consumption: < 0.1 mA

Working current consumption: < 500 mA

medicine Capacity: 6 mL

Average Nebulization Rate: $\geq 0.2\text{mL/min}$

Particle Size (MMAD): 3 to 5 μm

Aerosol output rate: 0.087 mL/min

Sound level: $\leq 21\text{ dB}$

Automatic off: No medicine / after continuous working for 20 minutes

Operating Temperature Range: 10°C to 40°C (50°F to 104°F)

Operating Humidity Range: 10 to 95% RH

Storage Temperature Range: -20°C to 50°C (-4°F to 122°F)

Standard package : main unit, nebulization unit, adult Mask, child mask, mouthpiece, pouch, USB cable, rubber protective cap

Product size: 6 * 5.5 * 13 cm

PARTS DESCRIPTION

(1) Main unit

(2) On / off button

(3) Jack for USB cable

(4) Adult mask

(5) Child mask

(6) USB cable

(7) Rubber protective cap

(8) Mouthpiece

(9) nebulizer release button

(10) lock of medicine container

(11) Nebulization unit

(12) Medicine container

(13) Mesh


(14) Aerosol output nozzle

(15) medicine inlet hole

(16) Vibrator (inside in the medicine inlet hole, not visible)

(17) Battery compartment

(18) Electrodes

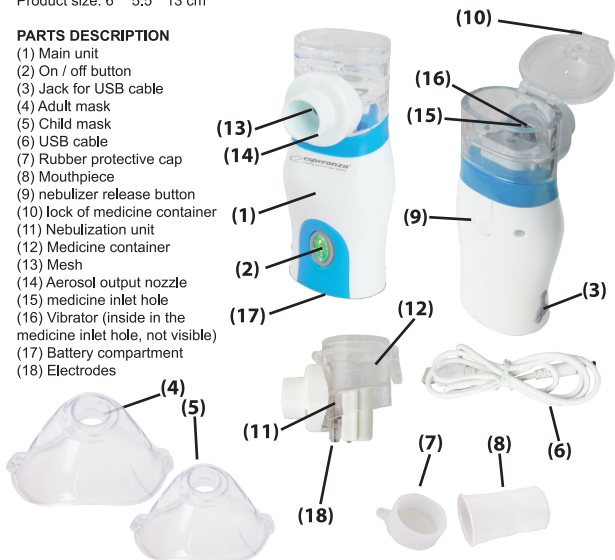HOW TO USE THE DEVICE

Connect to the power source

Select the power supply you prefer. You can connect the USB cable to the power source or use AA batteries.

USB cable: Connect one end of USB cable to the main unit and another end can be connect computer, laptop or adapter.

Alkaline batteries: The unit can operate with AA alkaline batteries. Open the battery compartment, follow the polarities marked inside the compartment to insert the new batteries and close the battery compartment.

How to use

- Press the medicine container release button and lift it.
- The mouthpiece, masks should be cleaned and disinfected before using the unit for the first time or after the unit has not been used for a long period of time, or daily during normal use. How to clean and disinfect, you can refer to the part " How to clean" and " How to disinfect".
- Pull up the lock to open the medicine container.
- Fill the medicine container with medicine and close the lock.
- Push the nebulization unit to the main unit till it clicks to the place.
- Select the power supply you prefer as described above.
- Put the mask or mouthpiece to the aerosol output nozzle.
- Press the "on/off" button and green light is on. The device start to work.

Note:

- If the device keep working for continuous 20 minutes, it will switch off automatically after 3 times blinks in red. Press the "on/off" button to switch on again.
- When the battery is low, the device will automatically switch off with 3 times blinks in red.
- When the medicine is finished, the device will automatically switch off with 3 times blinks in red.

HOW TO INHALE

Caution: Only use type and amount of medicine as prescribed and instructed by your doctor.

- Slightly tilt the unit to the position of mesh. In this position, the vibrator is immersed in the medicine and nebulization will start after power on. When the vibrator has been immersed in the medicine, the unit can be used under any angle.

Caution: In some positions (e.g. upright, nebulization may stop after a short while. In this case, briefly tilt the unit again in order to re-immersed the vibrator in the medicine.

- Place the mouthpiece in your mouth or place the mask over mouth and nose. Always be calm and relaxed when you inhale. Breathe in deeply so that the medicine can reach deep into the bronchial tubes. Hold your breath briefly, then breathe out slowly while removing the mouthpiece from your mouth. Do not breathe too quickly. Pause when you feel to have a rest. You may synchronize the nebulization with your inhalation/exhalation cycles by using the unit in manual nebulization mode.

Note:

- If excessive medicine collects on the mesh, nebulization may stop. In that case, turn off the power and remove the mask / month piece. Absorb the medicine by means of a lint-free cloth.
- Do NOT poke the mesh with a cotton swap or a pin; the mesh may be permanently damaged.
- After inhalation always switch off the nebulizer by pressing the on/off button.

HOW TO CLEAN THE UNIT

Note: if the device is not cleaned and disinfected correctly and frequently, micro-organisms may remain in the unit, causing risk of infection.

Cleaning of medicine container

1. Remove the mask/ mouthpiece from the main unit.
2. Remove the medicine container from the main unit:
3. Open the medicine container cap, discard the medicine residues.
4. Pour clean cold tap water to medicine container to remove medicine left in the vibrator holes. But do NOT immerse any other parts of the nebulization unit in water except the inside of medicine container.
5. Pour clean cold tap water again, close the container clock. Press the "on/off" button to start nebulization for 2 minutes. Pour away the excessive water. Let the container air dry on a clean, lint-free towel or paper.

Cleaning of nebulization unit, mesh and main unit

- Do not touch the mesh with a cotton swab, or with a pin.
6. Clean the main unit, mesh, main unit and electrodes by wiping them with a (moistened) gauze. Do not use volatile agents, such as benzene or thinner.

HOW TO DISINFECT THE DEVICE

Note: if the device is not cleaned and disinfected correctly and as frequently as indicated, microorganisms may remain in the unit, causing risk of infection.

- The mouthpiece, masks should be disinfected before using the unit for the first time or after the unit has not been used for a long period of time, or daily during normal use.
- Disinfection can be done by using alcohol. Soak the masks and mouthpiece in 70% ethyl alcohol for approx. 1 minute. Shake off excessive alcohol from the parts and let them air dry in a clean, lint-free towel or paper.
- Disinfection can also be done with vinegar. Use one part white vinegar diluted with three parts distilled water. Make sure mixed solution is enough to submerge mouthpiece and masks. Shake off excessive solution from the parts and let them air dry in a clean, lint-free towel or paper.

WARRANTY

1. This warranty is issued by EDC Poterek Sp. Jawna with address:
ul. Poznańska 129/133, 05-850 Ożarów Mazowiecki, Phone: 022-7213571 w.10 – called hereinafter Esperanza.
2. This warranty confers upon you specific rights and remedies against Esperanza. Please read this document carefully, by making a claim under this warranty, you will be deemed to have understood and accept its terms.
3. Esperanza warrants that any Esperanza product together with attached accessories is free from defects in materials and workmanship under conditions of its normal use.
4. If you make a claim under this warranty, Esperanza will, at its option, repair or replace any defective product into new one by using new or use factory-refurbished parts. In case of using factory-refurbished parts they will be equal to new products.
5. Defective product should be returned together with receipt, commercial invoice or any other proof of purchase to the place of purchase. It should include all accessories attached in original packing. Esperanza has right to reject warranty if any accessories will be missed.
6. Esperanza is committed to respond to the complaint within 14 calendar days of its notification in the point of sale. Esperanza will repair the defective product in the shortest possible time. If will be needed to import from abroad of any parts which are necessary to repair this time may be extended.
7. In case of replacing defective product onto new one, Esperanza has right to replace it with newer model with the same technical parameters if the old one model was replaced by new one during warranty period.
8. If you make claim under this warranty, you are deemed to have agreed that repair and (upon availability) replacement, as applicable, under the warranty services described herein, is your sole and exclusive remedy vis-a-vis Esperanza.

EXCLUSIONS:

This warranty does not extend to:

1. Any product not manufactured by Esperanza
2. Any product that has been damaged or rendered defective as a result of:
 - use it in other way than normal intended use
 - modification or internal interferes or changes of product
 - service or repair by anyone other than Esperanza service
 - improper store, transportation or packing
 - improper installation of product or software
 - mechanical, physical, chemical etc. damages
3. In no event shall Esperanza be liable for any damages, loss of business, profits, data or loss based on use Esperanza products with other devices.
4. Products with removed warranty and other labels and products spent or have hard use marks.

All products are covered by Esperanza's 2-year manufacturer's warranty while product is purchased by the consumer (the consumer is a person who buys product for the purpose of not professional or commercial activity). In the case of purchase for the company (for purposes relating to the activities of professional or commercial) warranty is 1 year from date of purchase. Esperanza products have all the necessary safety certificates required in the EU.

According to Directive 2002/96/EC, disposal of waste electrical and electronic equipment together with household waste is prohibited. Their components need to be recycled or disposed separately. Otherwise, they can cause the release of hazardous and noxious substances which pose a threat to the environment. The law requires each consumer to return unnecessary or waste electrical and electronic equipment to the producers, sales agents or waste treatment facilities. Detailed conditions are governed by national legislation. The above symbol placed on the product, in the user manual or on the packaging refers to this requirement. Thanks to proper segregation and compliance with the provisions relating to waste disposal each user contributes significantly to environmental protection.

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We, the manufacturer/importer (My, producent/importer):

DONGGUAN SIMZO ELECTRONIC TECHNOLOGY CO.,LTD

NO.6, Zhangzhou Road, Daojiao Town, Dongguan City, Guangdong Province, 523187, P.R.China

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Item No (numer/symbol): ECN005

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The products is compatible with the following norms/standards:

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Manufacturers/Authorised representative

 Jack Jin

MODEL IDENTITY DECLARATION

We, the manufacturer:

DONGGUAN SIMZO ELECTRONIC TECHNOLOGY CO.,LTD

NO.6, Zhangzhou Road, Daojiao Town, Dongguan City, Guangdong Province, 523187, P.R.China

hereby declare that following product:

Product description: MESH NEBULIZER

Model Number: NBM-2

Is technically identical with the following product imported from our company by EDC Poterek Sp. Jawna company:

Product description: MESH NEBULIZER

Model Number: ECN005

Trademark: ESPERANZA

Manufacturer authorized representative

Gray Jack Jin

EC Certificate Production Quality Assurance System: Certificate CN19/41053

The management system of

Dongguan SIMZO Electronic Technology Co., Ltd.

No.6, Zhangzhou Road, Daojiao Town, Dongguan City,
Guangdong Province, 523187, P.R. China

has been assessed and certified as meeting the requirements of

Directive 93/42/EEC

on medical devices, Annex V

For the following products

Infrared Thermometer used for measuring body temperature
(Model: HW-1, HW-2, HW-2S, HW-3, HW-4, HW-4S, HW-302, HW-303, HW-F1,
HW-F2, HW-F3, HW-F4, HW-F5, HW-F6, HW-F7);
**Portable (Ultrasonic) Nebulizer on Conscious Patients (NBM-1, NBM-2,
NBM-3, NBM-3B, NBM-4, NBM-4B, NBM-5, NBM-5B,
NBM-6, NBM-6B, NBM-7, NBM-7B, NBM-8, NBM-8B);**
Clinical Digital Thermometer
(TH-801, TH-802, TH-803, TH-804, TH-805, TH-806, TH-807, TH-808, TH-809)

Where the above scope includes class IIb or class III medical device(s), a valid EC Type Examination Certificate according to Annex III is a mandatory requirement for each device in addition to this certificate to place that device on the market.

This certificate is valid from 01 February 2021 until 24 May 2024
and remains valid subject to satisfactory surveillance audits.

Issue 3. Certified since 27 January 2015

and first certified by SGS Belgium NV since 16 December 2019

Certification is based on reports numbered CNDGGI 16426

Authorised by

Global Medical Devices Certification Manager

SGS Belgium NV, Notified Body 1639

SGS House Noorderlaan 87 2030 Antwerp Belgium
t +32 (0)3 545-48-48 f +32 (0)3 545-48-49 www.sgs.com

LFMD9009 - Certificate CE1639 Annex V_EN rev. 01

Page 1 of 1

This document is issued by the Company subject to its General Conditions of Certification Services, unless otherwise agreed, accessible at www.sgs.com/terms_and_conditions.htm. Attention is drawn to the limitations of liability, indemnification and jurisdictional issues established therein. The authenticity of this document may be verified at <https://www.sgs.com/en/verified-clients-and-products/verified-clients-directory>. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.